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
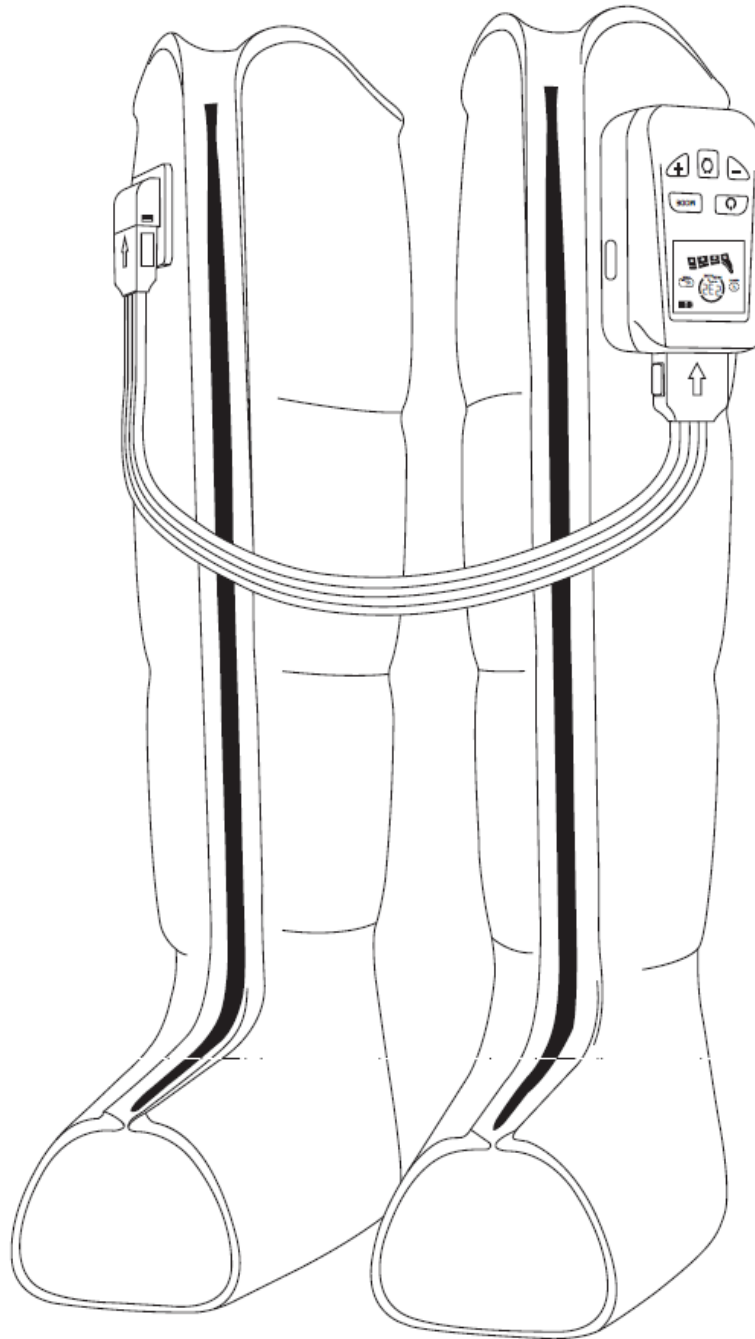
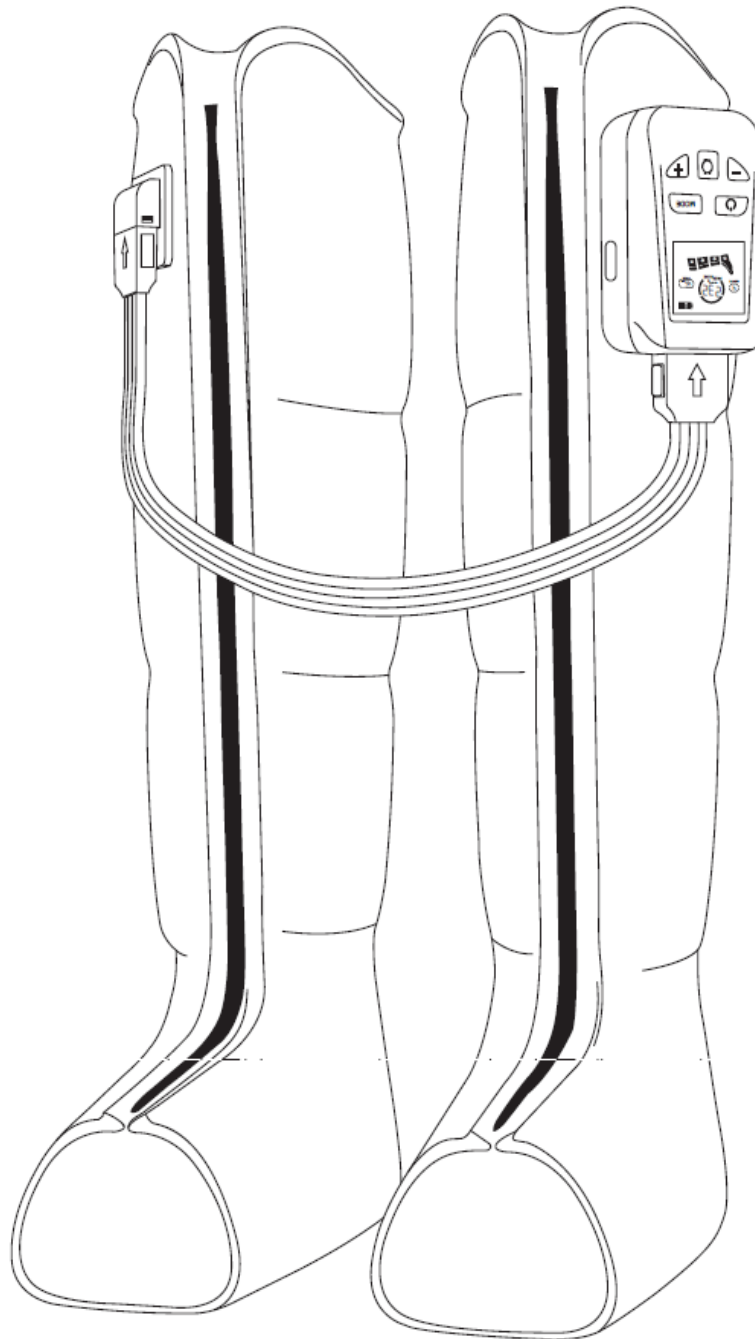
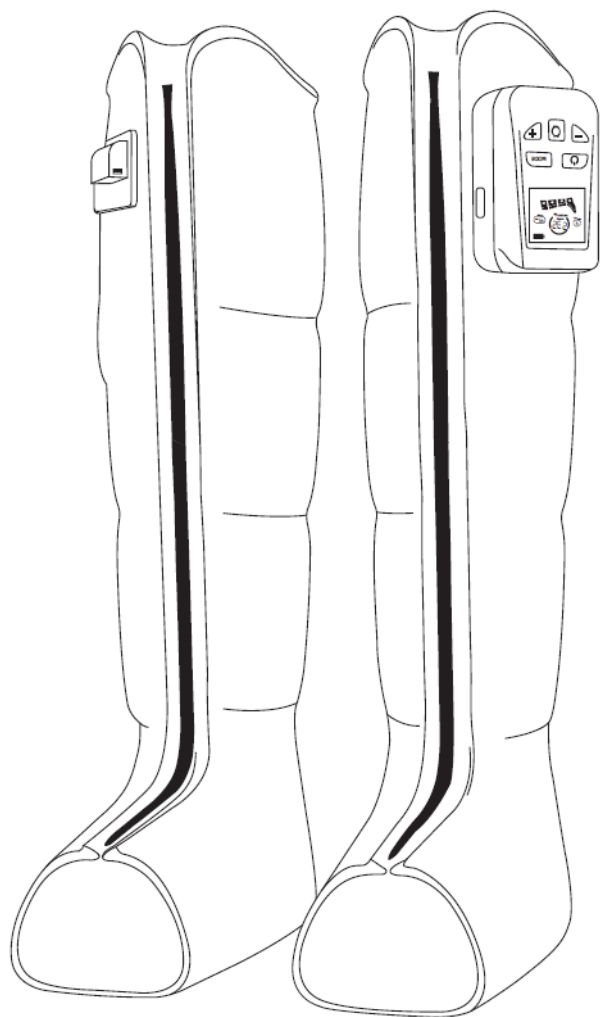
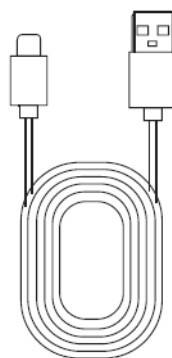
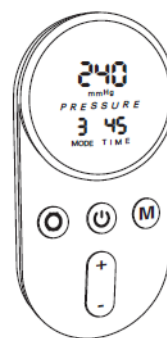
Remote control x 1 unit

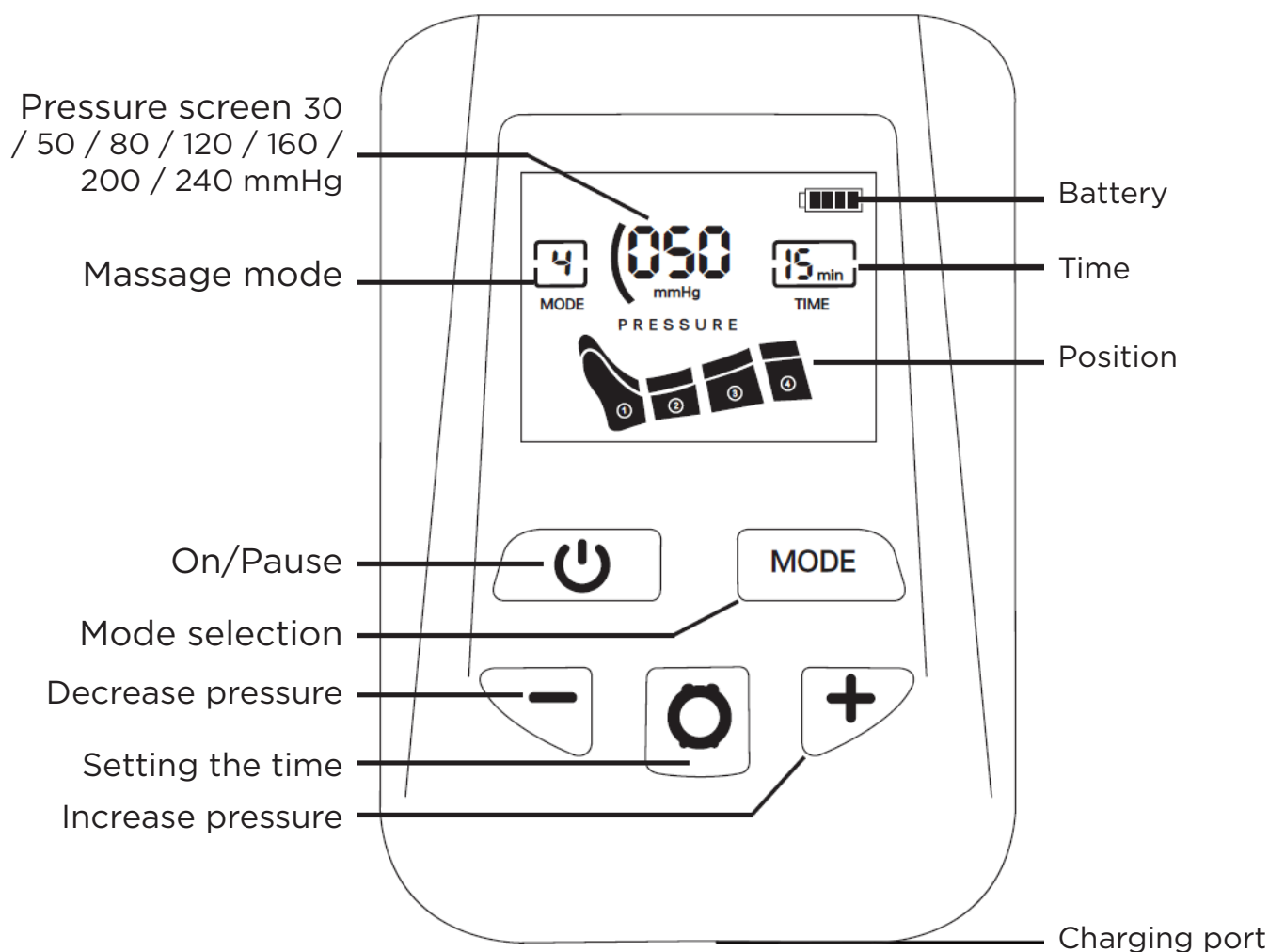Operating instructions


Step 1: Press the power button  on the LED panel or the button  on the remote control.

Step 2: Set the pressure by pressing the +/- button on the panel or remote control (30/50/80/120/160/200/240 mmHg)

Note: For first-time use, it is recommended to start with a low pressure, perhaps 80 mmHg, then you can adjust the pressure to the level that makes you feel most comfortable.


Step 3: Select the time by pressing the button .


Step 4: Select the mode by pressing the MODE button. The MODE mode will cycle through the following modes: (1, 2, 3, 4, 5).

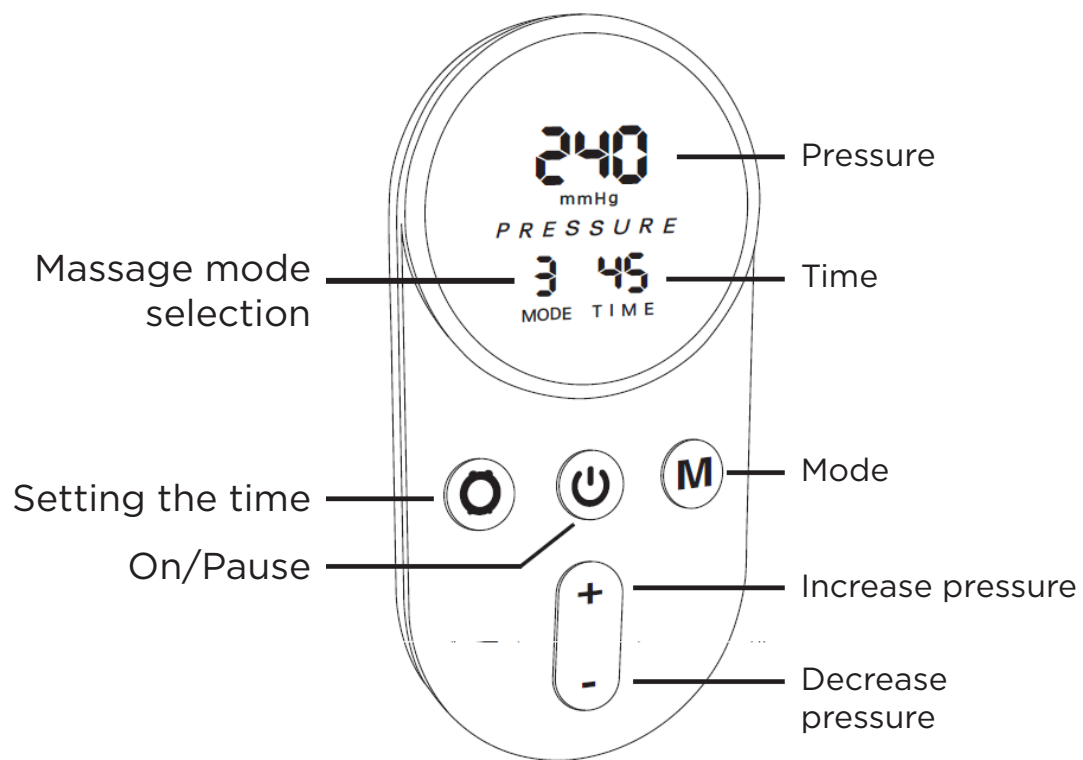Remote control (optional)

Illustration of the connection

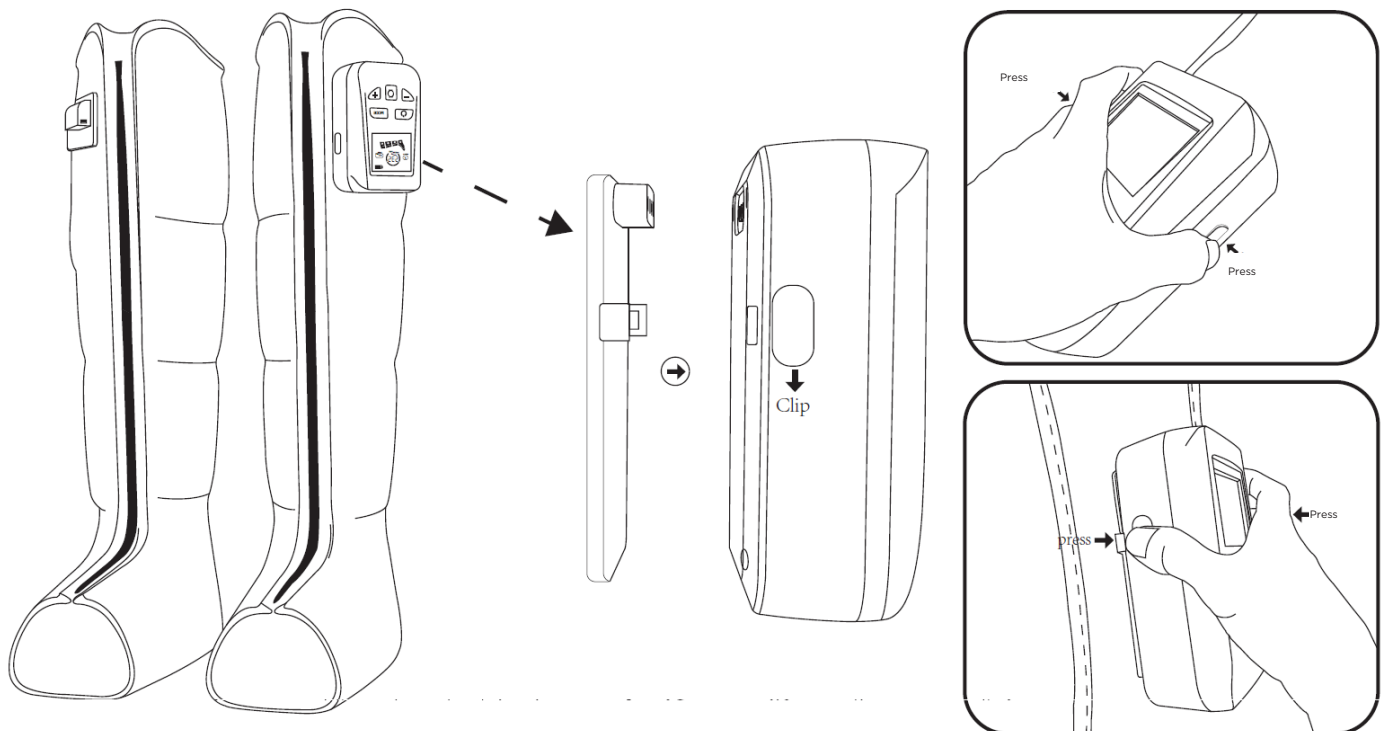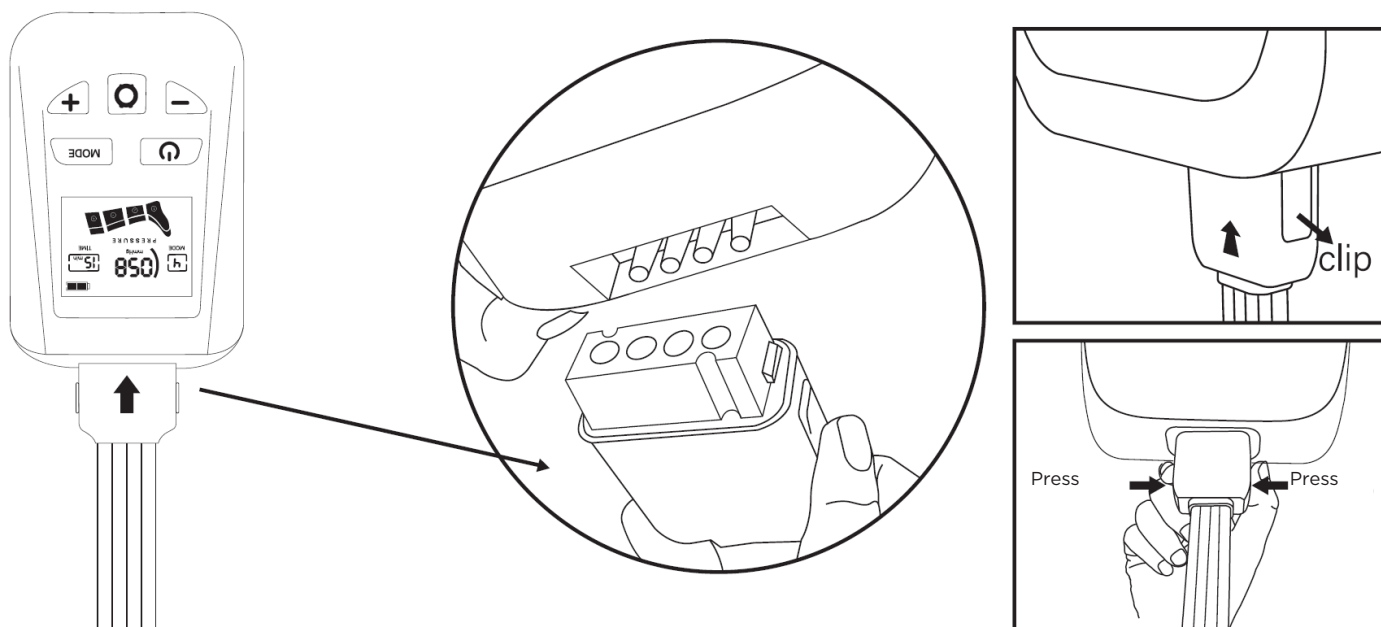After using the pump, when storing it, be careful to remove the hose cap by pressing the clips on the double sides of the hose cap to avoid breaking it.

Massage mode

Sequence: Automatic circular massage from bottom to top.

Circulation: Optional massage zone for individual massage.

Entire: Massage the entire leg gradually from bottom to top.

Wave massage

sequential compression

Double wave

All squeezing

Sequential squeezing + Double wave

Product introduction

The Air Compression Leg Massager is a smart massage device specifically designed to relax and restore leg muscles. It uses air pressure technology to simulate manual massage techniques, providing you with a complete and comfortable massage experience.

Equipped with multiple chambers, this massager covers the entire leg, including the thighs, calves, and ankles. As you use the massager, the chambers inflate and deflate sequentially, generating moderate pressure to promote blood circulation and relieve leg fatigue and muscle tension.

Battery

- Green indicator: 100%-65%
- Blue indicator: 65%-25%
- Red indicator: 25%-5%
- If the power level drops below 5%, the device will automatically shut off.
- CH will be displayed during charging. Please note that this product cannot be used while charging.

Product parameters

Pressure Level	30/50/80/120/160/200/240 mmHg
Noise	≤55 dB
Charging Port	DC
Battery	3200mAh/7,4V ; Lithium Battery
Charging Current	DC 12V-2A
Charging Time	2 hrs
Working Time	1,5-2 hrs
Product Weight	1,4 kg
Input Voltage	100-240V AC
Power	24W

Cleaning

Controller:

Clean the screen with a soft cloth.

Fabric:

This product cannot be washed. Use an alcohol sterilizer or wipe it with a clean, damp towel.

Storage requirements

Keep this product out of reach of children, away from moisture, heat, direct sunlight, and water. If not used for an extended period, it should be charged every two months.

Discard:

Dispose of discarded batteries according to local or environmental regulations.

Use environment

	Use environment	Storage and transportation
Temperature range	5°C ~ 40°C	- 10°C ~ + 60°C
Humidity range	5% ~ 93% RH	5 % ~ 90 % RH
Barometric range	86 kPa ~ 105 kPa	86 k Pa ~ 105 k Pa

Product failure analysis and treatment

Phenomenon	Possible causes	Solution
Product cannot start	Product battery depleted	Charging for a period of time before turning on
Product cannot start	Product battery depleted	Charging for a period of time before turning on
	If the timer is set, the machine will automatically turn off when the time is up.	Restart
Display Er 1 character	over voltage alarm	After automatic shutdown and voltage reduction, restart the machine
Low pressure display value	The product has air leakage	After-sales service can be requested within the warranty/useful life period

Product safety precautions

1. Charging with adapters that do not meet the requirements of this product is prohibited.
2. Unqualified personnel and those not part of the technical service are prohibited from disassembling and repairing this product.
3. If the power cord is damaged, to avoid hazards, it must be replaced by the manufacturer or its maintenance department and similar full-time personnel.
4. Do not use in humid areas or where the product has been ingested with water or other liquids.

5. After use, place the massager in a dry place, away from moisture and direct sunlight.
6. Do not place it in a flammable, explosive place with chemical solutions or damp.
7. Do not touch the product with sharp objects to avoid damaging the product parts.
8. If you feel uncomfortable or abnormal during use, discontinue use immediately.
9. Do not use in high-temperature environments.
10. Keep the product clean and use a dry towel during maintenance.
11. Before using this product, it is recommended that you consult a doctor or use it under the guidance of a professional:
 - Pregnant women and high-risk children should not use it.
 - Patients with mental health conditions
 - Patients suffering from skin infections, skin swelling, scald burns, etc.
 - High-risk patients who are seriously ill and unsuitable for massage
 - Injured people with weak or damaged bones
 - People with pacemakers
 - Persons who have other implanted metallic devices, such as prostheses
 - People with serious heart disease or abnormal blood pressure

Note: There is a possibility of personal injury and damage if used incorrectly.

Handle fragile items with care

Read the instructions carefully before using this product.

Moisture-proof and rainproof

Do not litter carelessly. Dispose of trash according to local regulations.