
GADNIC

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BALANZA INTELIGENTE CORPORAL

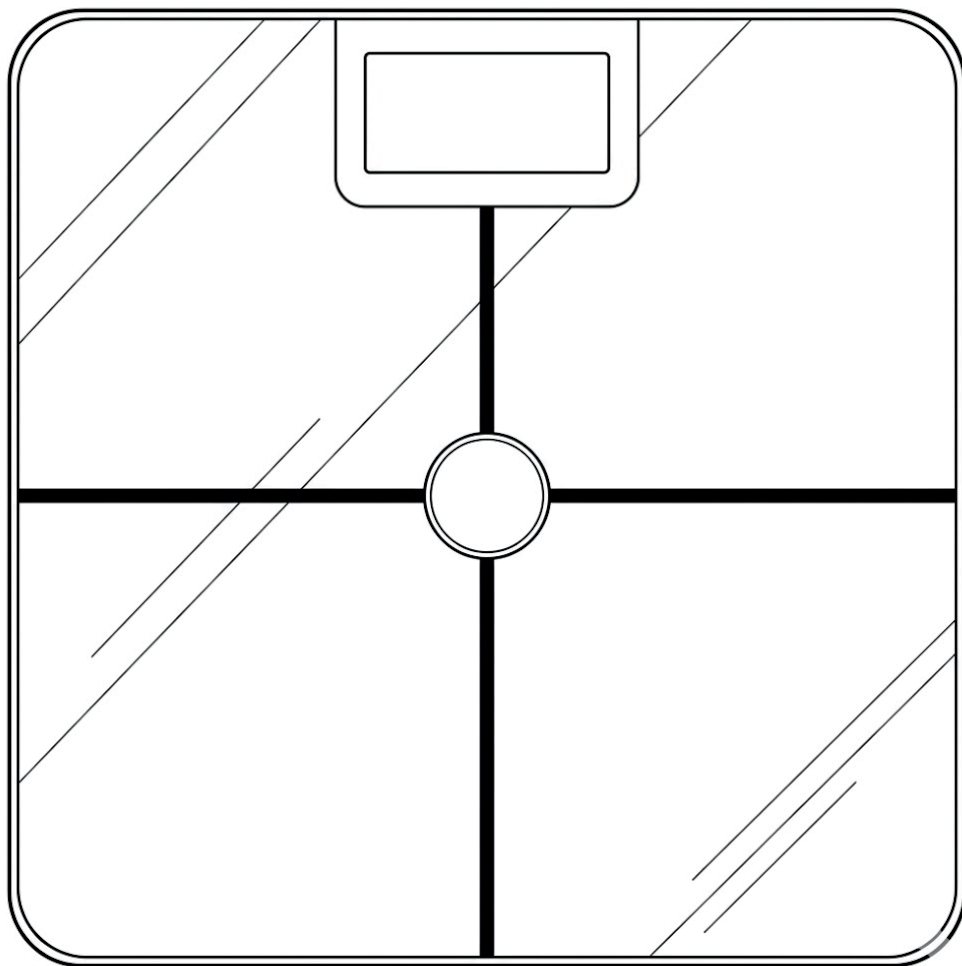
BALAN42X

INDICE

ESPAÑOL.....	03
ENGLISH.....	12
PORTUGUÊS.....	21

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BALANZA BALAN42X

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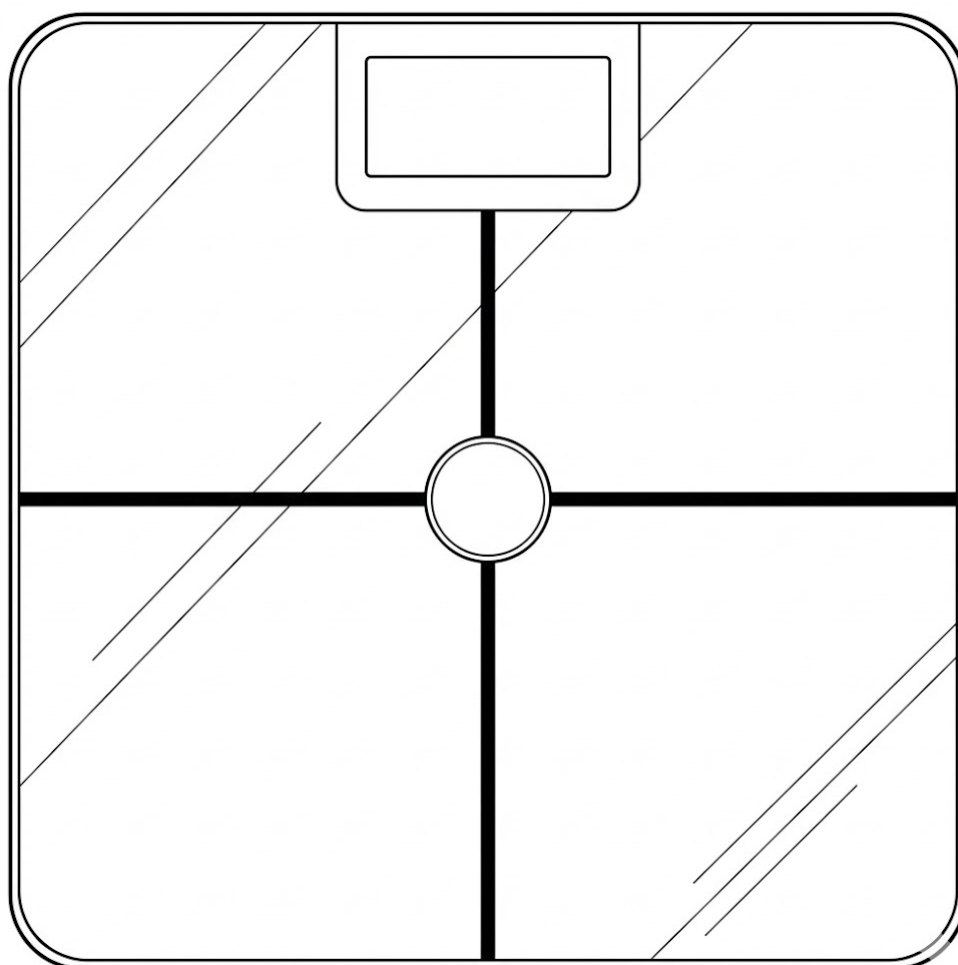
INSTRUCTIONS FOR USE

SCALE BALAN42X

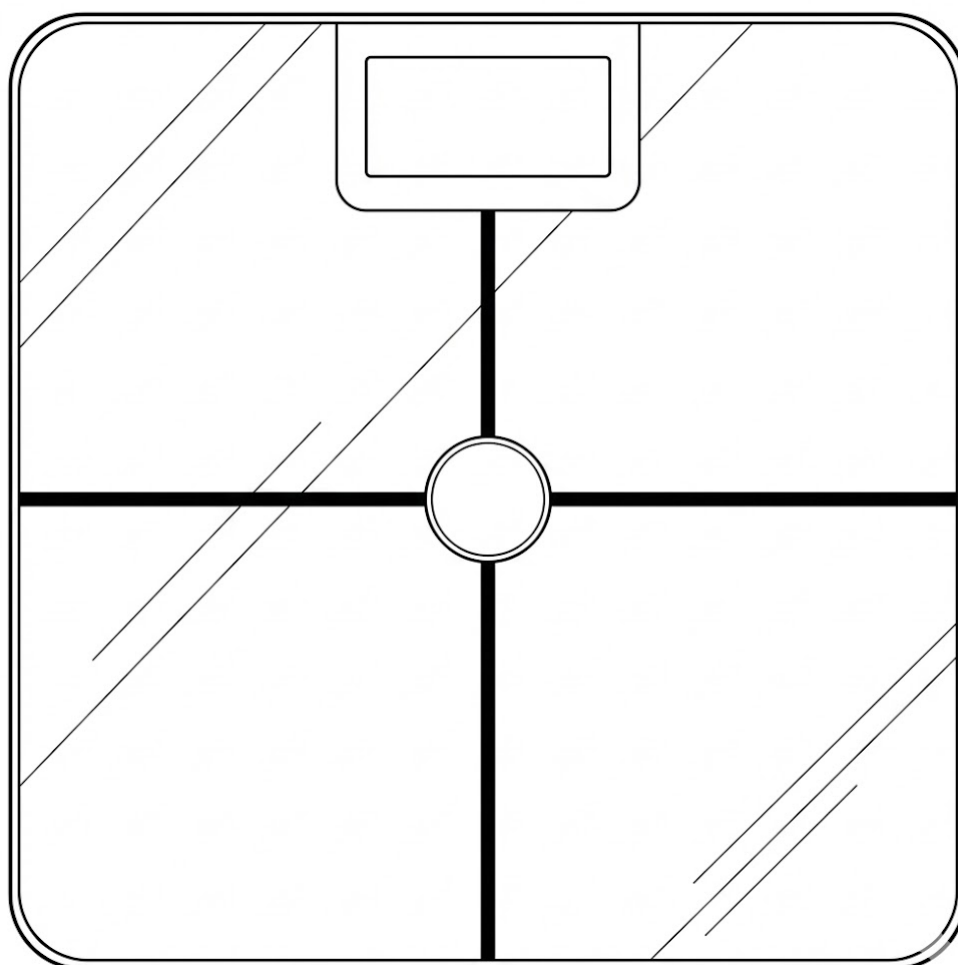PLEASE READ THE INSTRUCTIONS CAREFULLY BEFORE USE. KEEP THE MANUAL
FOR FUTURE REFERENCE

Product introduction

1. Use of body fat

To provide measurement of weight, body fat, protein, and other health data for people concerned about health, weight loss, and shaping.

2. Declaration

This device is not a medical device, and the parameter information and measurement analysis provided are for reference only. They are not used as a basis for diagnosis, prevention, treatment, or monitoring of diseases, and are not suitable for medical institutions such as clinics and hospitals.

The product illustrations, accessories, user interfaces, etc., in the manual are schematic diagrams for reference only. Due to product updates and improvements, there may be slight differences between the actual product and the schematic diagram. Please refer to the actual product for accuracy.

Power On

The battery-style body fat scale requires opening the battery cover and installing the battery in the correct positive and negative directions. After installation, the device will turn on automatically.

The charging-style body fat scale needs the button on the base (usually on the back near the screen) to be pressed slightly, and the electronic scale can be turned on by pressing it slightly.

Return to Zero

After turning on the electronic scale for the first time and placing it on the ground, the scale will show a weight of approximately 1KG, which is the weight of the electronic scale itself.

Método 1: Method 1: After the electronic scale shows the weight, it will show the reset

symbol “” for about 5 seconds and then show 0.00. At this time, the scale can be loaded normally.

Method 2: Step gently on the electronic scale up to a higher weight, and the scale will show 0.00, which can also be used for normal weighing.


If weighed directly without returning to zero, the weight shown on the electronic scale is your weight plus the weight of the electronic scale; Please ensure that the electronic scale shows 0.00 before using it.

To use the body fat scale product, it is necessary to use a smart APP to record and analyze your body data.

Search for “LEA ONE” in the app market/APP Store to download and install the application.

After the download is complete

1. To open the application for the first time, you need to authorize a location permission and Bluetooth permission. Please click “Allow”. If you skip this step, you need to open it in the application permission management of the phone settings application management. Only after enabling the corresponding permission can the mobile application bind to the body fat scale device normally.
2. Enter the login page
 - Login with account binding; Supports linked login with Apple ID, FACEBOOK account abroad, and WeChat account. Domestic users are recommended to use email registration login.
 - It should be noted that the data obtained from user measurements will be linked to your login account and stored on the cloud server. Please use the most common login method. If you need to register and log in using an email, go to the upper right corner of the login interface. Click register to log in using your email.

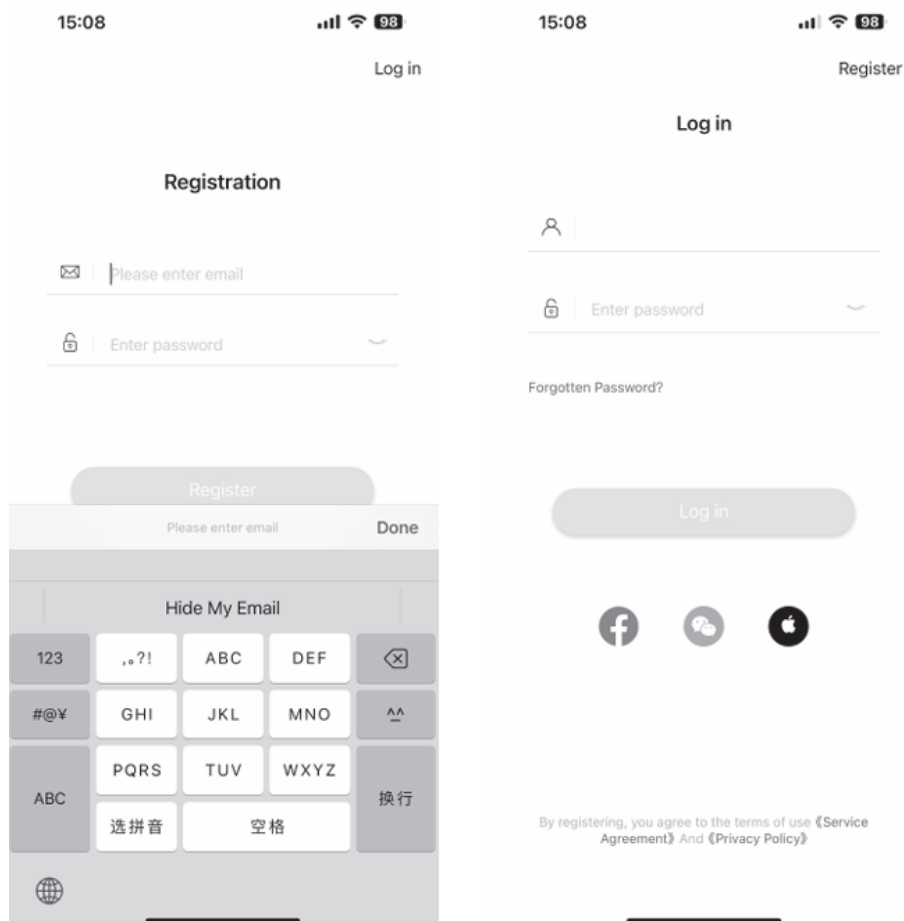Device Binding

1. Before starting device binding, it is necessary to turn on the body fat scale. Battery-powered devices require the AAA battery to be inserted into the battery compartment.
2. Place the electronic scale on a hard, flat surface, and enter the Add Device page in the software; Click 'Add Device' and step lightly on the body fat scale to change the reading on the display screen. You can see the device number and icon in the "Devices nearby" tab, click the icon to complete the binding.
3. After completing the device name, click 'Start using' to enter the measurement page.

✓ Matched successfully


Start using


Measurement Instructions

- To ensure measurement accuracy, try to measure at the same time and under the same conditions, and avoid measuring in the following situations:
 - After intense exercise.
 - After taking a shower or sauna, eating excessively, having a cold or fever, drinking a lot of water or alcohol.
1. Please place the product on a hard, flat floor. Do not place it on carpets or other soft surfaces as it may cause inaccurate measurement data.
 2. Before initial use, please step lightly on the scale surface and reset it to zero to eliminate its own weight.

1

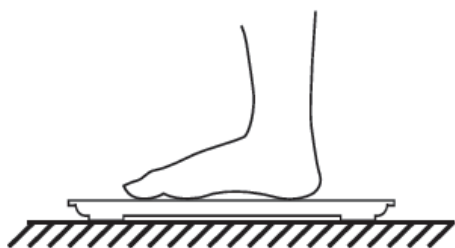Place the scale on a hard and flat floor

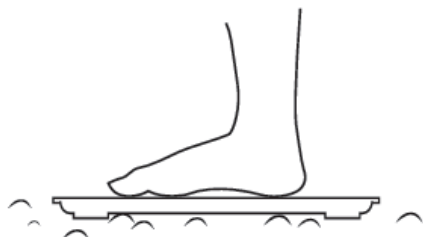Do not place on carpets or soft floors

2

Scale measurement barefoot, waiting for the weight at the end of the scale. When the number remains fixed, it indicates that the weighing is complete.

3


Second Measurement

1. Under the Equipment tab > Click on the device icon that has already been linked.
2. Enter the previous data page > Shows “Connecting...” at the top of the page
3. Repeat steps 1-2-3 in the measurement instructions.
4. You can update the data immediately.

Other Functions


1. Clicking on a specific data item can only view the reasonable range of the data; And various data explanations.
2. Click on the “---” in the top right corner to achieve data sharing, data comparison, and data sharing.
3. Click on the avatar in the top left corner to add family members, up to a maximum of 7 family members can be added.

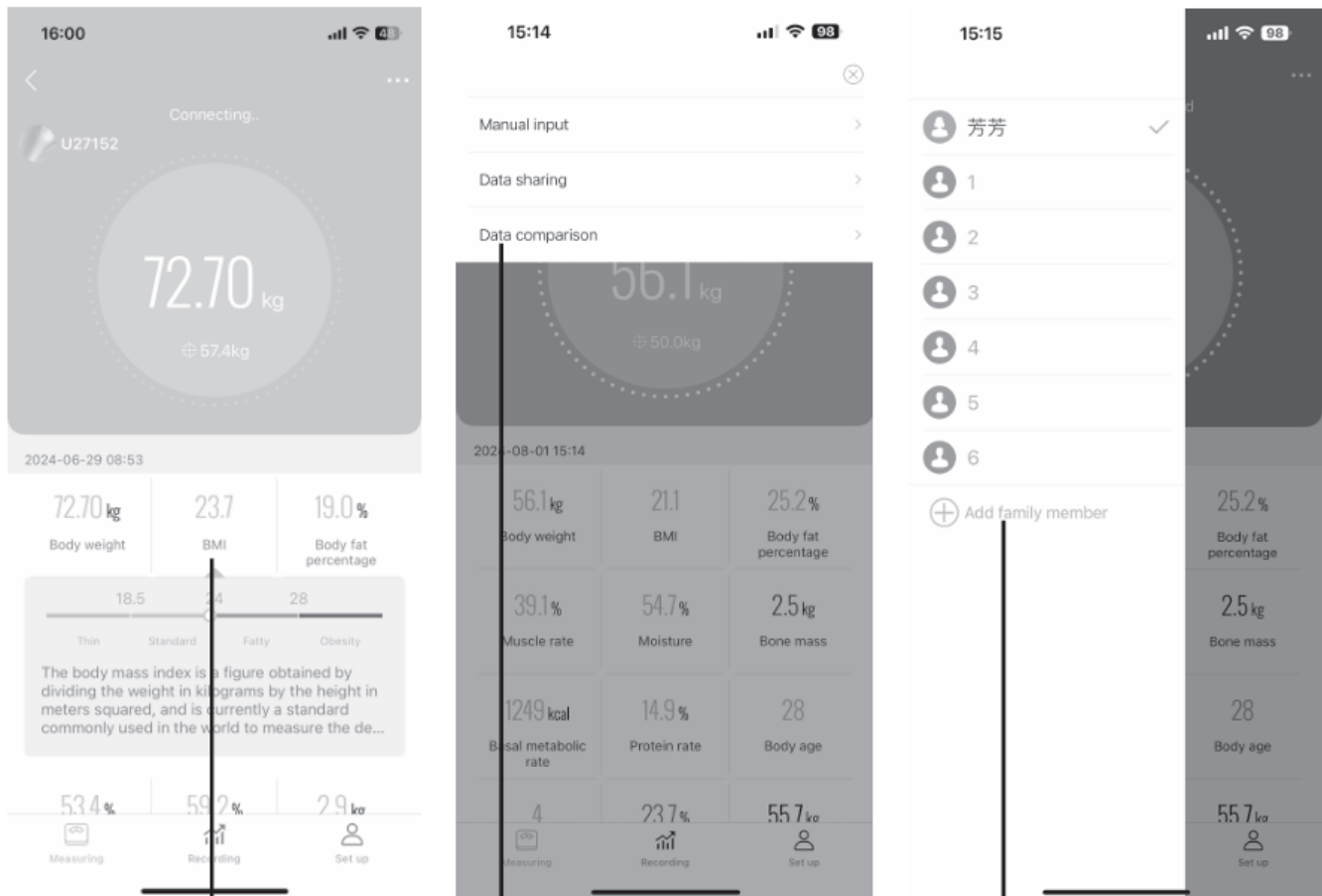After entering family members, manual switching is required before starting measurement. Family members, followed by measurement.

The data comparison function allows for a comprehensive comparison of changes in various data at different times.

Click on the record to view the curve chart of various data, which can provide feedback on the long-term trend changes of your body.

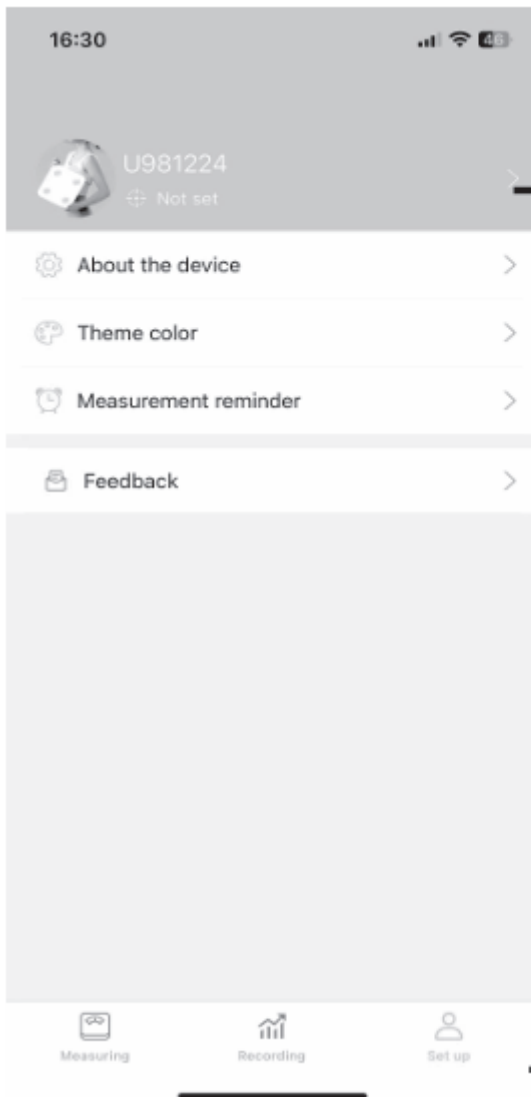Click on the avatar to modify your own basic data: Height, age, gender, nickname.

Settings Tab


Matters that need attention

- To ensure data accuracy, please complete the measurement maintaining your feet dry and clean, wearing light clothing, barefoot, and maintaining the skin of your feet in close contact with the metal electrode pads.
- Do not move your feet or leave the scale surface at will during measurement.
- Please measure at the same time and under the same conditions every day.
- Please try to avoid using it after intense exercise, excessive alcohol consumption, excessive intake of water or food, sauna bathing, and wait for the body to recover before taking measurements to avoid affecting the measurement results.
- Do not submerge this product in water or other liquids.
- Do not weigh yourself on the scale when your feet or the surface are wet to avoid slipping.

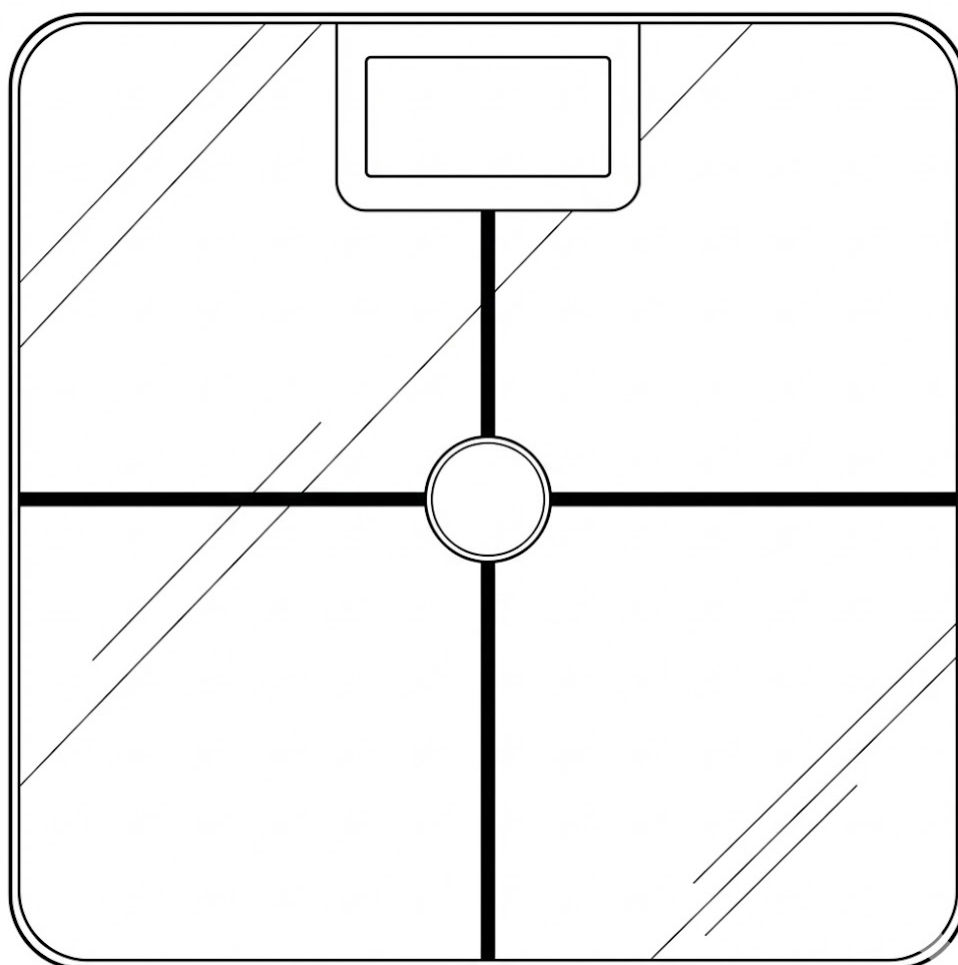- Do not use this product in high-temperature environments above 40C, below 0C, or near strong magnetic fields.
- Keep away from humidity, water sources, fire sources, and environments with a lot of dust and corrosivity.
- Keep away from electric heaters, microwave ovens, ovens, water heaters, stoves, candles, and other places that may generate high temperatures to avoid affecting battery performance and causing equipment failure.
- Do not crash, squeeze, drop, throw, or scratch this product. Strong impacts may cause partial breakage of the scale surface.
- If the product is damaged, please stop using it immediately and consult after-sales service promptly. Do not touch or attempt to remove broken parts.
- Please operate this device in accordance with the recommended instructions for use in this manual.

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