
# GADNIC

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### BALANZA DIGITAL ULTRA FINA

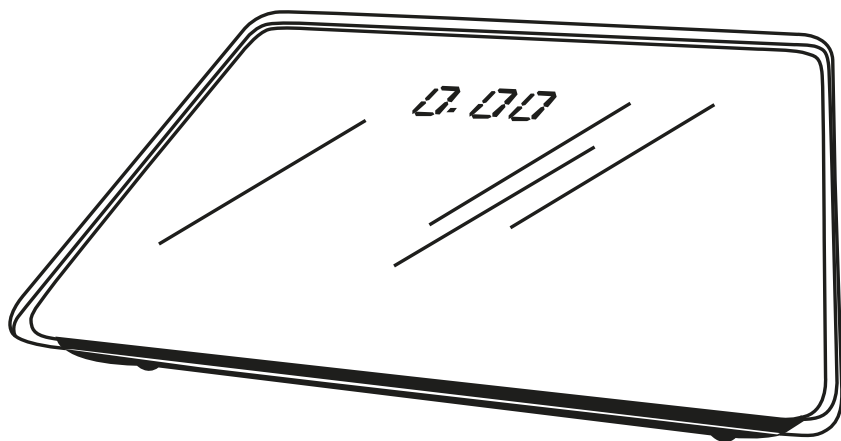**BALANZ16 / BALANZ17**

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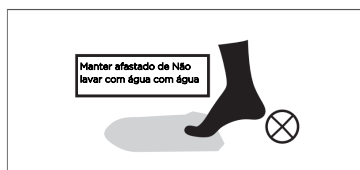
## Specifications

- ★ Starting weight: 5kg
- ★ Capacity: 180kg
- ★ Graduation: 50g
- ★ Display: 4-digit LCD
- ★ Power supply: 2x1.5V (AAA) - Batteries not included
- ★ Operating temperature: 5-35 degrees Celsius
- ★ Weight unit: kg
- ★ Automatic switch-off: 10 seconds
- ★ Low voltage indication: 'Lo'
- ★ Overweight indication: 'oL'.

## Instructions

- ★ Place the scale on a flat surface and step on it when the display shows 0.00 kg. shows 0.00 kg.
- ★ Stand firmly in the centre of the tray. When the weight value is stable, the weight value is locked.
- ★ If you leave the tray and step on it again after it is locked, it can be weighed again. weighed again.
- ★ If you don't step off the tray after locking it, the weight value on the tray will remain unchanged for 10 seconds. The weight value on the pan remains unchanged for 10 seconds before automatically switching off.
- ★ If the scale shows 'Lo', the batteries are low and should be replaced.

## Attention

- ★ Do not leave heavy objects on the scale for long periods of time when it is not in use. when the scale is not in use.
- ★ When weighing, try to stand in the centre of the scale.
- ★ Stand firmly on the scale.
- ★ Do not hit the scale with force,
- ★ Do not use the scale in the presence of water.
- ★ Keep the scale and your feet dry during weighing.

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