*Care*  
by GADNIC

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**EXFPIE01**

## INDICE

ESPAÑOL.....	03
ENGLISH.....	10
PORTUGUÊS.....	17

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## REMOVEDOR DE CALLOS RECAR- GABLE EXFPIE01

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# INSTRUCTIONS FOR USE

## RECHARGEABLE CALLUS REMO- VER EXFPIE01


PLEASE READ THE INSTRUCTIONS CAREFULLY BEFORE USE. KEEP THE MANUAL FOR FUTURE REFERENCE

# Product Introduction

Congratulations on your purchase! This Rechargeable Callus Remover is an innovative device created to smooth dry, rough and calloused skin on the feet. unique rollers grinding spins 360 degrees, at an amazing 40 times a second, 300 degrees wide angle to gently and effectively buff away hard, jagged skin on the heels. This Callous Remover works faster and is safer than metal scrapers, with no blades to cut the skin, and easier to use than other methods, that require constant rubbing to achieve results, it is battery-operated and designed to fit comfortably in your hands, smooth, silky feet has never been easier.

This Callous Remover is recommended for one person only, for hygienic reasons, unique rollers grinding should not be shared with others.

## General descriptions

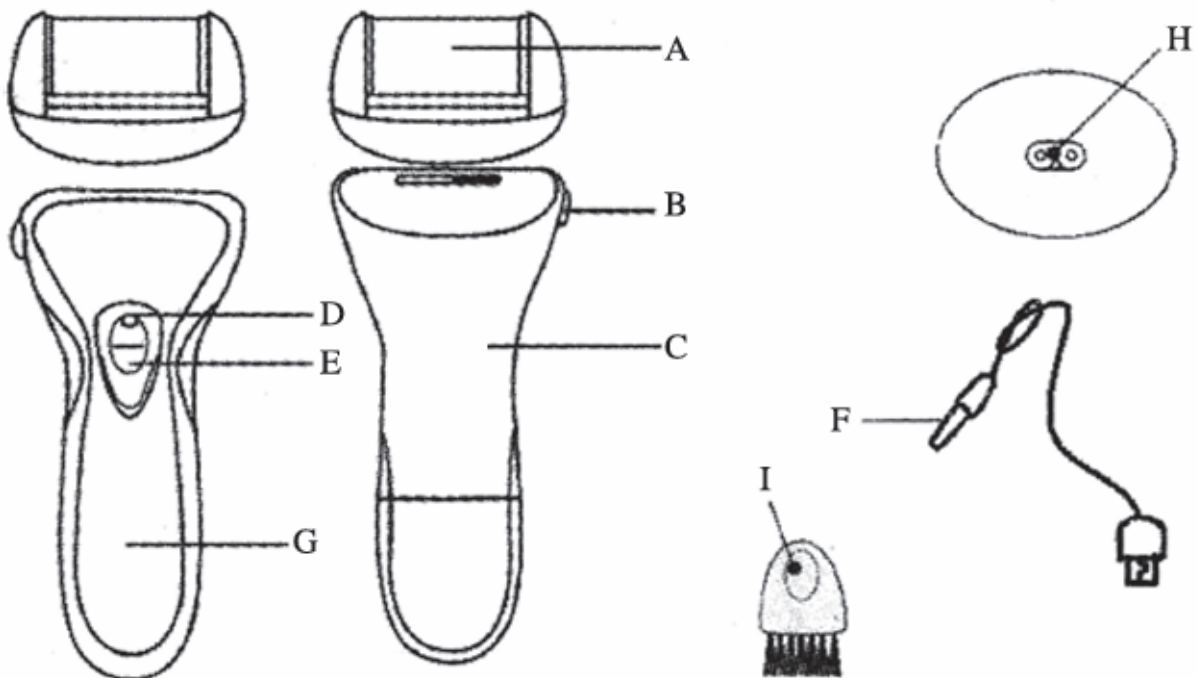- A. Unique Rollers/Grinding Head
- B. Release Button
- C. Anti-Skid Handle
- D. Safety Switch Lock
- E. ON/OFF Switch

- F. Power Adaptor
- G. Indicate Light
- H. Charging Socket
- I. Cleaning Brush


# Matters needing attention

1. The Callous Remover is only for use on the feet and heels. Do not use this device anywhere else on the body.
2. Use the Callous Remover only for purposes described in this instruction book.
3. Each head is recommended for one person only. For hygienic reason, the head should not be shared with others.
4. Do not use if you have diabetes or poor blood circulation.
5. If you have any skin conditions or other medical issues, please consult with your physician prior to use.
6. Do not use on skin for more than 2-3 seconds at a time.
7. Do not use on skin that is red, irritated, inflamed, infected or blistered.
8. Do not use on warts, on sores or open wounds.
9. Discontinue use immediately if you experience any pain or irritation.
10. Overuse may cause skin irritation.
11. Seek medical attention if skin becomes severely irritated or bleeds.
12. Household use only.
13. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning the use of the appliance by a person responsible for their safety. **Warning/Important**
14. Do not place or store your unit in direct sunlight. Please store your unit in a cool dry place at a temperature between 0-40°C. Keep the appliance dry, when it is not in use.
15. Do not turn your unit on if the rollers grinding head is not attached.
16. Do not drop or insert anything into any opening of the unit.
17. Do not use your unit outdoors.
18. Children should be supervised to ensure that they do not play with this appliance.
19. Keep your unit out of reach of children as it is not intended for use by children.
20. This appliance is not to be used on, by or near children or invalids.
21. To prevent accidents and injury to you or damage to the unit - when your Callous Remover is in use, keep it away from the hair on your scalp, your eyebrows, eyelashes as well as clothes, brushes, wires, cords, shoelaces, etc.

22. Do not continually rub the unit on the same area or apply strong pressure.
23. Keep the protective cover on your unit when it is not in use.
24. Do not disassemble this device.
25. Products regularly charging time for 4 hours.

## Charging Instruction

1. Make sure the on/off switch is in the off position (Fig 1).
2. Insert the USB power cord into the appliance charging socket.
3. Insert the USB charging cord into the power plug and then connect with power supply or computer USB interface to charge (Fig 2).
4. The charging light is lit up, make is charging.
5. The electricity charge indicator will be always on.
6. Before using this appliance, or appliance is not used for long time, should be charge 4 hours.
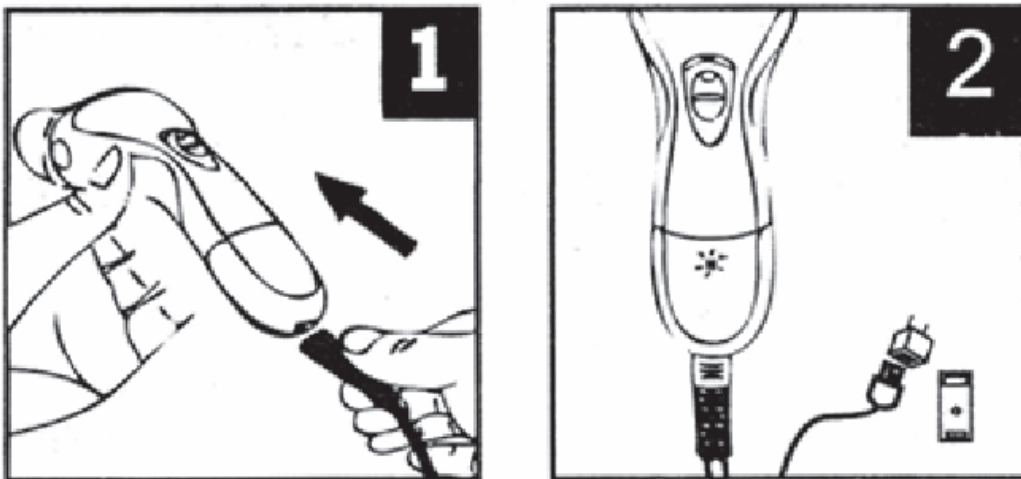
## How to use the Callus Remover

This Callous Remover is designed for smoothing coarse skin on the feet and heels. Do not use this device anywhere else on the body.


1. Wash and clean the skin and dry fully. Note: Do not use on broken or bleeding skin.
2. Make sure the Rollers Grinding Head is placed into the unit. See section of Changing the Rollers Grinding Head on how to remove and replace head.
3. Remove the protective cover from the unit and it is ready to Use.
4. Switch on (fig. 3).

5. Gently rotate (back and forth or side to side) (fig. 4) the callous/corn/hard skin area for 2-3 seconds. Do not press hard-gently glide the unit over desire area. The unit will stop if it is press too hard on the skin. Stop and check if you have achieved the softness you desire. If not, go over the area for 2-3 more seconds and check again. Continue to use this method until the level of smoothness is obtained.
  - Warning: stop using if skin becomes sore or Inflamed.
  - Warning: Do not use on skin for more than 2-3 seconds at a time.
6. Turn the unit off (fig 5).
7. Rinse the skin or use a wet towel to get rid of the excess dead particles
8. To clean this Callous Remover, please see the Cleaning section.
9. Keep the protective cover on your unit when it is not in use.

## How to use Nail Polisher

1. Repeat the above steps, install the polisher head into the top fo unit.
2. First use the buffer head to buff the nail first (fig. 6). Do not use it on the same place for long time when in buffing.
3. After buffering, clean the surface of nail.
4. Use the polisher head to make the nail shinning (fig. 7).

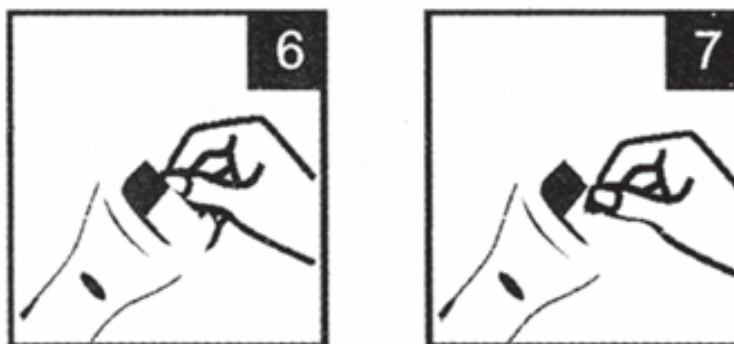
# Changing the Rollers Grinding Head

Change the Rollers Grinding Head as needed for best result. the Rollers Grinding Head last for a few sessions depending on frequency of use and size of area being treated. Replace as needed. Make sure the unit is turned off before replacing the Rollers Grinding Head.

1. One hand hold the head, another hand press in the release button, on the side of the unit, the head will tack out (fig. 8 & 9).
2. Place a new Rollers Grinding Head into the unit until you hear a “click’ sound. Be sure to insert the new head into the unit, it has been firmly installed.


Note: Do not turn your unit On if the Rollers Grinding head is not attached.
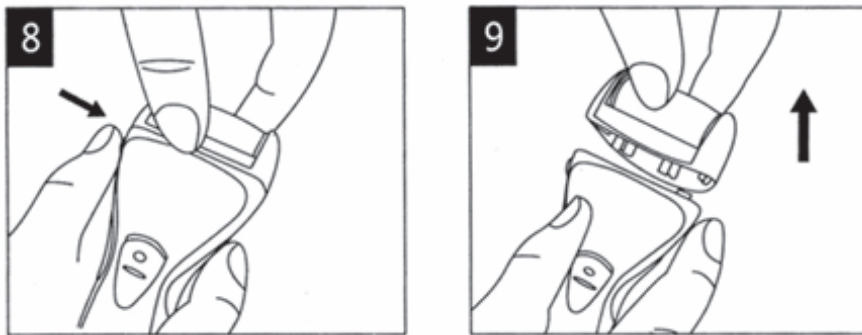
## How to Cleaning

Note: Clean the appliance after each use to ensure better performance. In order to clean this Callous Remover the safest way, follow the instructions below:

1. Remove the protective cover from the head.
2. Clean the Callous Remover with the cleaning brush. Do not switch the appliance on while brushing (fig. 10).
3. Clean the exterior of the unit with a slightly damp cloth to remove debris. Do not place the unit under the tap or submerge it in water.
4. Ensure the unit is turned off.
5. Wipe the body of the unit with a clean dry cloth.
6. Replace the protective cover.

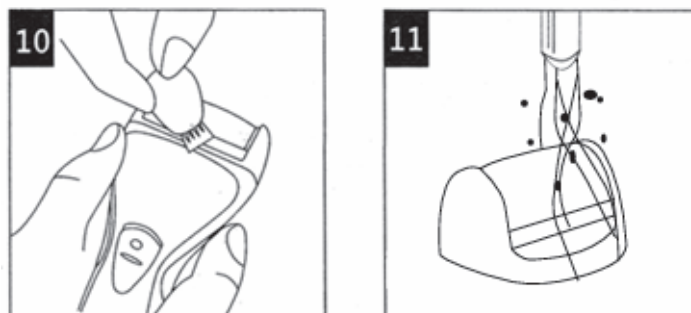**Warning: The callus remover is not washable, please keep it dry.**

## **Storage**

- When the appliance is not in use, put it on a dry place that children can not touch.
- Do not place the appliance in an area with excessive heat or moisture.

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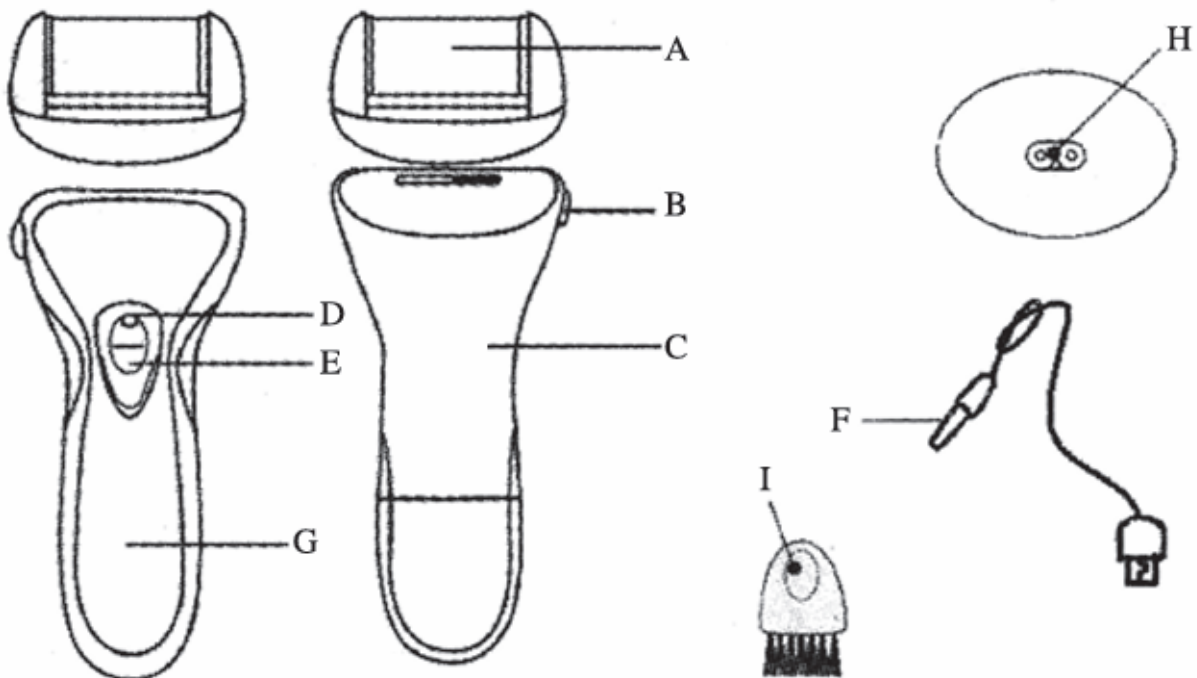
E. Interruptor Ligar/Desligar (ON/OFF)

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