GADNIC

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MASAJ001

INDICE

ESPAÑOL.....	03
PORTUGUÊS	09
ENGLISH.....	15

PRECAUTION

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PRECAUTION

Thank you for purchasing our Easy Health Kneading Massage Cushion. Before using the Easy Health Kneading Massage Cushion, please read the User Manual carefully to learn the correct methods of use, especially the safety notes.

FUNCTION AND SPECIALTY

1. 3 MAIN FUNCTIONS

Shiatsu, Kneading and Vibration modes.

2. 3 TYPES OF MASSAGE MODE

Upper back massage, lower back massage and full back massage.

3. 3 TYPES OF KNEADING MESSAGE PROGRAM:

Massage up and down, massage the tips and knead the neck.

4. INDEPENDENT CONTROL


For seat vibration massage.

5. FULL FUNCTION REMOTE CONTROL.

MAINTENANCE


- 1.** Turn off and unplug the power cable before cleaning the unit.
- 2.** Use a soft, damp cloth to clean the pad.
- 3.** Do not use thinners, volatile oils or bleaching powders that cause the fabric to change color or become damaged.
- 4.** Do not use lubricant on this product.
- 5.** If the cushion is not in use, keep it in a dry and less dusty place.

PRODUCT

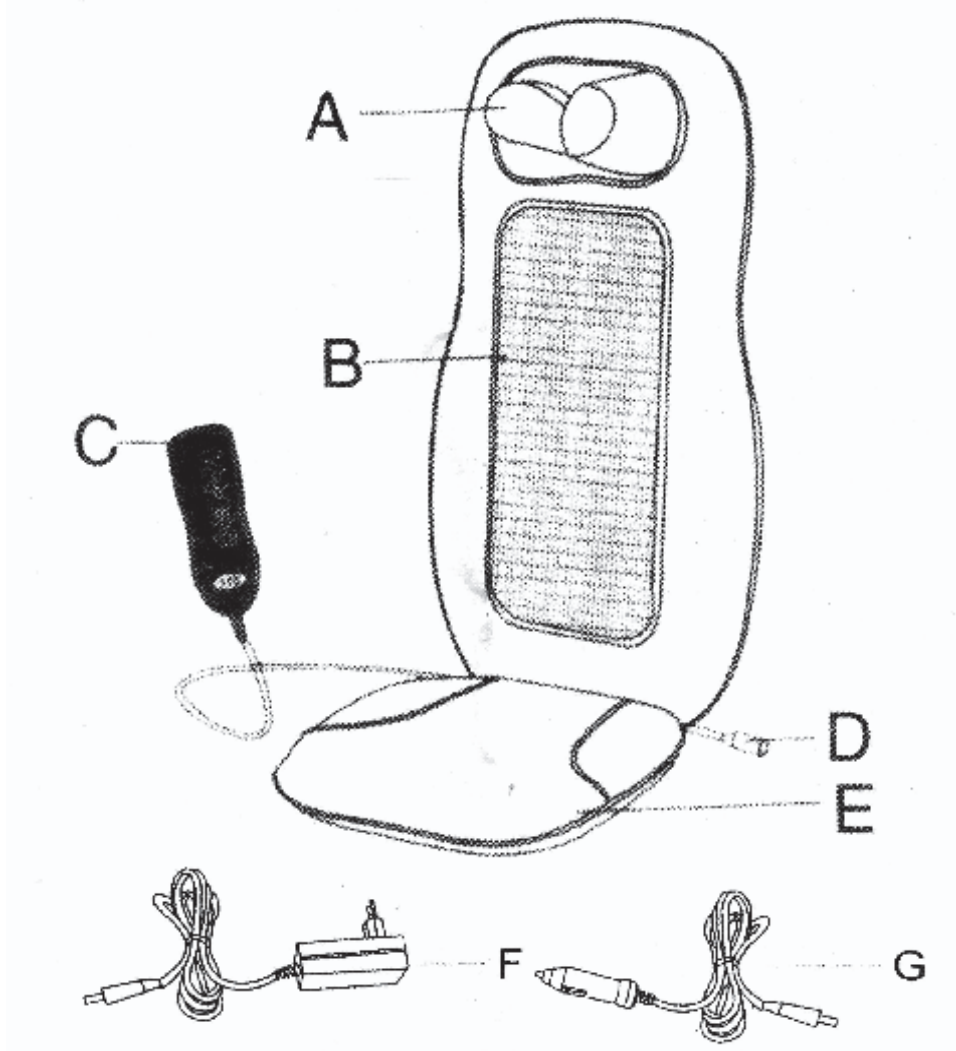- A. Neck massage cushion.
- B. Back massage cushion.
- C. Remote control.
- D. Power connector for the massager.
- E. Vibration seat
- F. Power adapter.
- G. Car adapter

REMOTE CONTROL

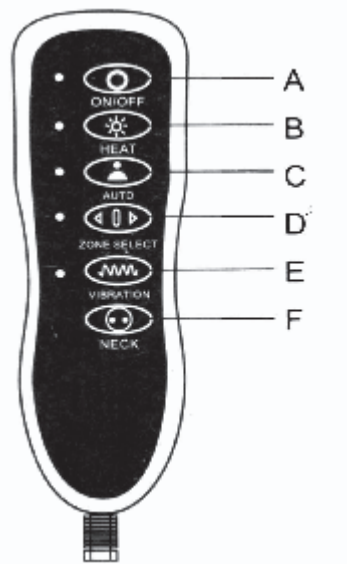- A. On / Off
- B. Heating
- C. Selecting the massage area
- D. Vibration
- E. Neck kneading massage

METHODS OF USE

- 1.** Connect the power adapter to the charging cable and plug the power adapter into the mains socket.
- 2.** Place the massage cushion and accessory on the chair.
- 3.** Press the “On/Off” button and the cushion will start operating in the Auto program.
- 4.** Press the “HEAT” button to make the infrared massage head glow and heat up; press the button again and the infrared and heat will be switched off.
- 5.** Press the “Select Massage Area” button and the “Auto” program will be changed to “Manual”.
- 6.** Press the AUTO button when the manual program is in progress; it will switch to Auto mode.
- 7.** Press the “Vibration” button. The seat vibrator will turn on; press it again and it will turn off.
- 8.** Press the “Neck Knead” button and the massage head will start working. Press the button again and it will turn off.

PRODUCT MAINTENANCE

- 1-** If you hear the sound “GiGi” coming from the cushion, this is the normal sound of the motor.
- 2-** If the sound is very intense during the massage, it may be that the cushion has been used too much. Turn off the cushion for at least 15 minutes. After that, you can turn it on again.
- 3-** If you can’t turn on the massage cushion, check to see if there is a warm sensation on the cover of the unit. If there is, the problem is caused by the overheating cut-off system. We suggest that you turn it off and wait an hour before turning it on again.

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