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NECK AND SHOULDER MASSAGER ·
MASSAJADOR CERVICAL DE
PESCOÇO

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

NECK AND SHOULDER MASSAGER

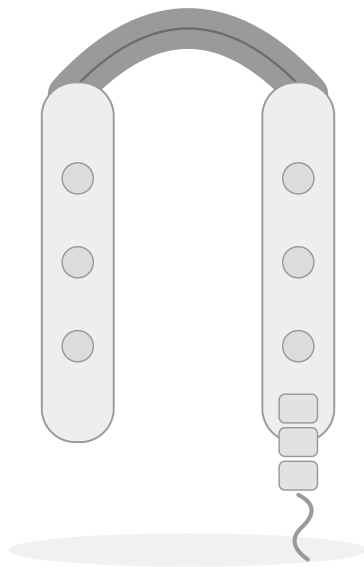Illustrative image

READ THE INSTRUCTIONS CAREFULLY BEFORE USE. KEEP THE MANUAL FOR FUTURE REFERENCE.

PARTS DIAGRAM

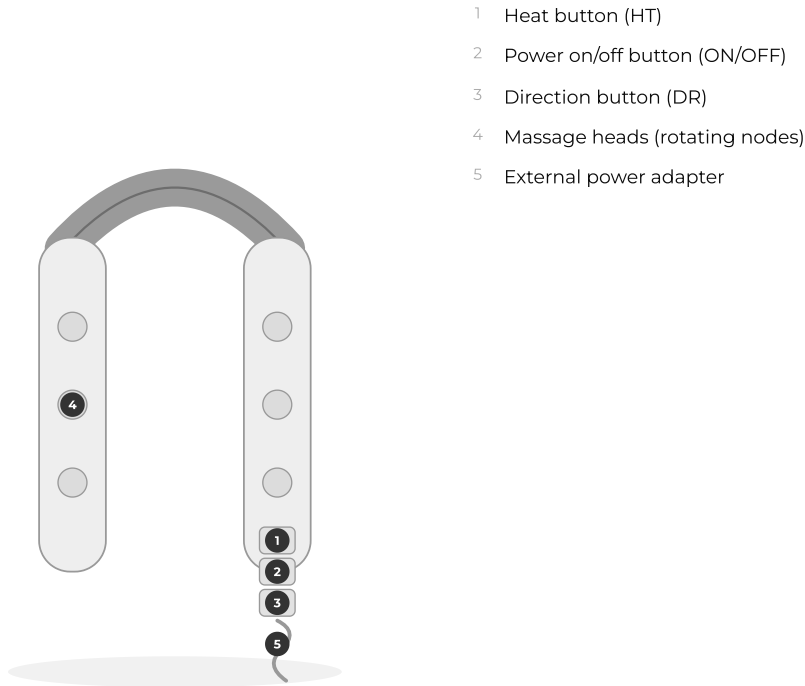SPECIFICATIONS

PRODUCT NAME	Kneading (shiatsu) neck massager
PRODUCT WEIGHT	1.1 kg
OPERATING TIME	15 minutes with automatic shutoff
POWER SUPPLY	Via external transformer/adaptor; check voltage compatibility before use
RECOMMENDED USE	Home and vehicle (car adapter included)

BENEFITS OF REGULAR MASSAGE

IMPROVES SLEEP QUALITY	Helps the body relax, eliminates fatigue, and relieves pressure.
IMPROVES BODY FUNCTION	Promotes blood circulation, aids lymphatic detoxification, and relaxes muscles.
SUPPORTS GASTROINTESTINAL FUNCTION	Aids digestion and relieves constipation and other digestive discomfort.

INSTRUCTIONS FOR USE

1. POSITIONING THE MASSAGER

Place the massager around your neck and shoulders, with the massage heads resting on the desired areas.

2. TURNING IT ON

Press the ON/OFF button to turn on the massager. It will shut off automatically after 15 minutes of continuous use.

3. ADJUSTING THE HEAT

Press the HT button to turn on the heating function; press it again to turn it off.

4. ADJUSTING THE MASSAGE DIRECTION

Press the DR button to change the rotation direction of the heads and focus the massage on specific areas of the body.

5. USING IT IN A VEHICLE

Connect the car adapter to a 12V power source to use the massager inside your vehicle.

DANGER

- After using the product or before cleaning it, unplug it.
 - If the product falls into water, do not touch the wet product; unplug it first.
 - Do not use the product while showering or bathing.
 - Do not place the product in water or any other liquid.
-

WARNING

- When no one is using it, unplug the adapter; if not in use, keep it unplugged.
 - To avoid overheating, electric shock, or injury, do not use the product on pillows or blankets.
 - Keep it away from children.
 - Use it according to the instructions in this manual.
-

NOTE

- Use only the transformer/adapter recommended by the manufacturer; confirm voltage compatibility before use.
 - Make sure the product and the transformer are properly connected.
 - Do not pull the connecting cable by force; a poor connection may cause the product to stop working.
 - Do not use it in the bathroom or in humid environments.
 - Do not place the product on hard surfaces while using it, as it could be damaged.
 - The product generates heat on its surface; people who are insensitive to heat should pay special attention.
-

TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSES AND SOLUTIONS
The massager won't start	Make sure the connector is properly seated in the switch. Confirm that the switch is working and that the button has been pressed correctly.
It stops working during use	Check whether the connector has come loose or there is a poor connection. You may have pressed the control button without noticing, or the preset 15-minute time may have elapsed.

If the problem persists, turn the product off and unplug it. Do not disassemble or attempt to repair it yourself; ask the seller or a service technician for help.

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