

## STORAGE / CARE

Regular care will extend the life of your massager also ensuring safe operation at all times.

- 1.Recommended storage: dry, free of dust.
- 2.Keep liquids and open reaway from the Massager
- 3.Clean the equipment using a dry, soft cloth.
- 4.Always unplug the Massager when not in use.

## CONFORMITYDECLARATION

Corresponds to the European guidelines 2006/95/EC.

## WARRANTY CONDITIONS

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt.

The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period.

The warranty claim does not apply if it is established that a defect has occurred e.g. through external in uence or as a result of repair or alteration not carried out by the manufacturer or an authorized dealer.


The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product.

If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer.

# GADNIC

## MANUAL DE USUARIO

### MASAJEADOR ELÉCTRICO MASAJ052

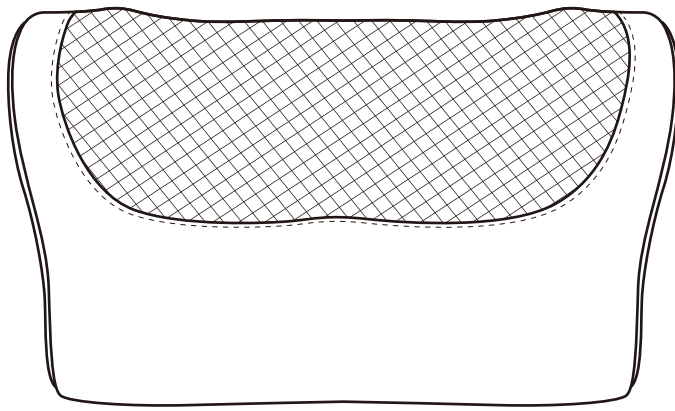

**MASAJ052**

## CONGRATULATIONS!

With the purchase of this massage device you have shown awareness of your health.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal Massager.

## TABLE OF CONTENTS

1. Safety instructions
2. Inside the package
3. Start up
4. Technical data
5. Equipment / functions
6. additional features
7. Control panel
8. Troubleshooting
9. Storage / care
10. Conformity declaration
11. Warranty conditions

## SAFETY INSTRUCTIONS



Please carefully read the following instructions to before using this tness equipment to en-sure its proper function and optimal e ect. Please retain these instructions for use!

message device complies with the recognized principles of technology and current safety regulars-

- This equipment is maintenance free. Any necessary repairs must be performed by an authorized pro-fessional. Improper use and unauthorized repairs are prohibited for safety reasons and will void the warranty.
- Never touch the power plug with wet hands.
- Please avoid the equipment coming into contact with water, high temperatures and direct sunlight.
- Be sure to unplug the equipment from the power point during extended periods of non-use to eliminate the risk of a short circuit.
- Do not use damaged plugs, cords or loose sockets. Damaged plug or cords must be replaced by the manufacturer, service representative or quailed personnel. In the event of malfunction unplug the equipment immediately. We assume no liability for damages resulting from wrongful or improper use.
- In order to prevent excessive stress on muscles and nerves, the recommended massage time of 30 minutes should not be exceeded.
- Never use pointy or sharp objects.
- Any massage – also massage by hand – has to be avoided during pregnancy and/or in case of one or more complaints in the area of the massage: head diseases recent injuries, thrombotic illnesses, in am-mations and swelling of any kind as well as cancer. Consultation of your physician is recommended prior to any treatment of ictions and ailments.
- Do not use the device while driving.
- If you depend on electrical aids such as a pace maker, please consult your physician prior to usage in any case.

## INSIDE THE PACKAGE

- Massager
- Adapter

## START UP

Please connect the adapter with the device and subsequently plug the power plug into the power socket. Now turn on the device by the control panel.

## TECHNICAL DATA

Measurements device:	40X25X9cm
Weight:	1.4kg
Power source:	
USA	input: 120V~50~60HZ 700mA / output: 12.0 VDC 2500mA
CANADA	input: 120V~50~60HZ 700mA / output: 12.0 VDC 1600mA
AUSTRALIA	input: 100-240V~50~60HZ 1000mA / output: 12.0 VDC 1600mA
EUROPA	input: 200-240V~50HZ / output: 12.0 VDC 1600mA
Rated power:	24 Watt

Certificates:



## EQUIPMENT / FUNCTIONS

The Massager's integrated deep warmth intensifies the massage. It permeates the massage heads and quickly becomes noticeable. During the warm up period, the circular movements of the massage heads relax tense muscles. The rotating movement of the massage heads loosen up the tissues.

Thanks to its ergonomic form the can give an effective massage to nearly all body parts-be it neck, shoulder, back, legs, feet or stomach. Thanks to the optimal placing of the massage heads with their relaxing circular motion, a soothing comforting massage is guaranteed every time.


Regular massage with the stimulates the functions of the nerves and activates the metabo-lism. Blood circulation is improved, even the most stubborn muscular tension is gently massaged away. In certain cases improved blood circulation and higher metabolism can even prevent blood clots and malfunctions of the internal organs.

Anyone who has had the opportunity to enjoy a Shi-atsu-Massage knows how beneficial it is. Regardless of the purpose of the massage, be it as a preventive measure against possible ailments or simply to relax and improve the quality of your life, a good, effective massage is an unforgettable experience.

The central philosophy of the Shiatsu Kneading is that the human being is a combination of body and spirit.This massage technique combines ancient traditional knowledge with modern treatments. The Shiatsu-Massage relaxes and regenerates mind and body giving you a new lease of life; relieving stress and at the same time giving you a feeling of invigor-ating well-being.

## CONTROL PANEL

METHOD 1: USE POWER ADAPTER

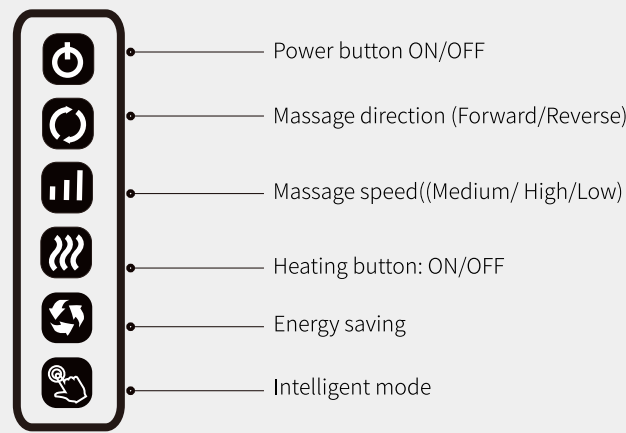

Step 1: Make familiar with the control buttons below.

Power button ON/OFF

Massage direction (Forward/Reverse)

Massage speed((Medium/ High/Low)

Heating button: ON/OFF

Note: the unit will be shut-off automatically in 15 minutes.

Stop using the product, please switch off the power button.

## TROUBLESHOOTING

Error

The message device cannot operate eventhough it has been connected to the powersupply.

Adjustment

- Check the connection of the adapter.

- The device could be overheated after prolonged use. It will automatically cool down to resume its normal condition. Please turn o the power and allow the device to cool o for 30 to 50 minutes.

# SERVITECH

SERVICIO TÉCNICO OFICIAL GADNIC  
[www.servitech.com.ar](http://www.servitech.com.ar)

☎ 11 6260 1114 (sólo texto)  
[serviciotecnico@servitech.com.ar](mailto:serviciotecnico@servitech.com.ar)