GADNIC

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MASAJ154

⚠ Advertencia

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model:	BC-8816/8817/8819	Appliance: 
Rated voltage:	12V	Rated time: 15min
Rated power:	24W	Net weight: about 2.0kg
Executive standard:	GB4706.1-2005 GB4706.10-2008	

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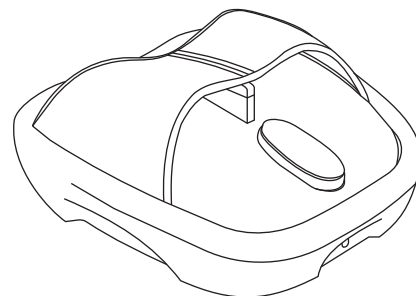
⚠ Warnings

To avoid burning, power short circuit or other hazards to users, be sure to read the following precautions.

1. When the product is not in use, please be sure to unplug the power plug.
2. Please do not stand on the product. Attention must be paid to people who are not sensitive to heat when using the product.
3. Please use the product only for the purposes mentioned in the.
4. Please keep the power cord away from heat source, do not use the power cord to pull the product, and do not let nails or wet objects into the product.
5. Don't use the product in the bathroom or other damp places, which may cause short circuit or fault of the product in the bathroom or other damp places, which may cause short circuit or fault of the product. Do not drop water on the product.
6. Do not use the product on heating body, such as electric blanket.
7. If the product is damaged, stop using it immediately to prevent accidental damage. It must be replaced at the designated maintenance point of the company before use.
8. Do not use damaged power cord or plug. If the power cord or plug is damaged, It must be replaced with a dedicated cord or a dedicated component purchased from its manufacturer or maintenance department.
9. Make sure to turn off the switch before unplugging the power plug. Do not use the product through the transformer, which may cause machine failure or short circuit.
10. In order to avoid short circuit or other faults, please do not disassemble any parts of the product by yourself.
11. Do not sleep while using the product. Do not use the product after drinking.
12. Do not allow children to play with the product as a toy. Children and the mentally handicapped should use the product under the supervision of others.
13. Do not unplug the power plug when the machine is running.

⚠ Warm prompt

★ Before starting the device to enjoy the massage, put the feet in the first place, put the foot forward as far as possible to massage, adjust the position and then press the open key.


★ Due to the different degree of personal massage force, you can wear thick socks or wrap dry towels at the soles of your feet to cushion the strength.

★ Applicable size: up to 43 yards.


If you have any of the following conditions or are taking any medication, consult your doctor before using it.

- It is forbidden for pregnant women, lactating women, infants and those with osteoporosis.
- Patients with thrombosis, aneurysms, dilated veins, or other conditions should not use massage products on doctor's orders.
- People with spinal curvature, heart disease and abnormal physical conditions.
- If you feel discomfort or pain when using this product, please stop using it immediately.

The panel control

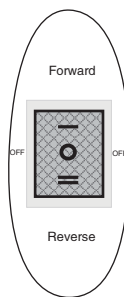(BC-8816)

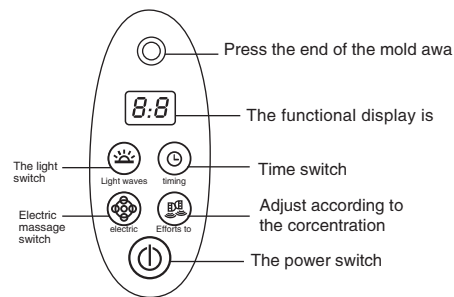(BC-8817/8819)

- 1.) ⏻ power switch: press this button once and the pedicure machine will enter the standby state, and the digital display will be 15 minutes(countdown).
- 2.) ⚙ Electric massage: electric foot massage on, off and reverse switch, the default off display: PO, automatic mode: P1, reverse display:P2
- 3.) ☀ Warm feet: Switch foot warming mode: default HO-off, H1-low, H2-medium, H3-high.
- 4.) ⚙ Power mode: switch ten modes: Default FO-off (F1-E9, Low-Strong
- 5.) ⚙ B8.8 function display: display various functions.
- 6.) ⌚ Timing function: switch the timing time, and the screen displays the set time, the unit is minutes BC-8816).

Ship switch forward; Massage motor began to work, toward the back press; Reverse work according to the mini-motor. Stop working in the middle.

Troubleshooting

I normal operation condition

- Please confirm whether the power adapter is plugged in and whether the power supply is smooth and the connection with the host is successful.

I The machine is suspended during use.

- If the operation stops after 15 minutes, the operation stops after the end of automatic massage mode.
- If multiple electrical products are used at the same time, the power will be too large, and the product safety device will automatically stop operating. In this case, the machine can be shut down and restarted. Do not stand while massaging

Cleaning and maintenance

1. Before cleaning the massager, remove the power plug first.
2. Do not use benzene, diluents or other corrosive solvents to clean your massager.
3. Do not use a wet cloth to clean the power switch of the massager, do not use alcohol, dry cotton cloth must be used.
4. If you do not use it for a long time, it is recommended to put away the power adapter and store the massager in a clean and dry environment.
5. Do not store the massager in a place with high temperature, open fire or direct sunlight.
6. Massage device should be handled with care, avoid falling, pressure, collision or strong impact.

To the user

Dear users:

Hello! Thank you for choosing our products. In order to ensure the correct use of the pedicure machine health care Products and give full play to its health care massage function, please read the instructions carefully when using. If you do not have any questions, please contact with our after-sale service department directly.

Thank you for your.

Cooperation!

Technical parameters

model:	BC-8816/8817/8819	Appliance: ☐
Rated voltage:	12V	Rated time: 15min
Rated power:	24W	Net weight: about 2.0kg
Executive standard:	GB4706.1-2005 GB4706.10-2008	

Warranty provisions:

This product is guaranteed for one year from the date of purchase (from the date of issuance of the sales voucher). All according to the provisions of the normal use of the fault, after our company confirmed, the company will be within the warranty period of the product for free maintenance, then the relevant freight borne by the customer.

Warranty card

Product model: _____ Order No.: _____
User name: _____ Purchase date: _____
Address: _____ Telephone: _____
Cause of failure: _____


Packing accessories pedicure machine, 1 manual, 1 power adapter

Warnings

Massage health • a
healthy life

Air Foot massage
Full of
pleasure
Massage pedicure
machine

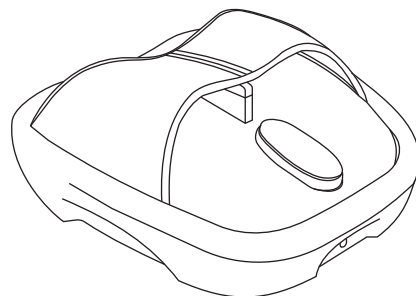Information in this manual is subject to change
without prior notice.

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