GADNIC

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MASAJEADOR ABDOMINAL PORTÁTIL

MASAJ196

Índice

Español	5
Português	12
English	19

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Mini masajeador TENS

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User Manual

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1. Product functions

1. Promote blood circulation and relieve body soreness and fatigue from working for a long time.
2. Helps the body regularly eliminate toxins from the body.
3. Pulse current is applied to muscle deep massage to achieve the effect of weight loss.

2. Product characteristics

1. High-tech body care products,LCD screen display, small current work, safe and reliable.
 2. Stylish appearance, simple, compact, easy to operate.
 3. Intelligent design, 8 modes and 19 intensities for you to choose from, single channel dual output.
 4. 8 modes and 19 intensities can be switched at will, without worrying about the intensity being too high or too low.
- Comfortably enjoy a variety of massage modes: massage, training, kneading, squeezing, beating, patting, shaping, activation, a total of 8 modes.

3. Specification


Device Name	Mini masajeador TENS
Model	BISD-MS9125
Power Supply	3.7V, 110mAh
Battery	Polymer Rechargeable Battery
Battery Duration	About 150min
Main engine weight	About 18g
Fully Charged	Approx. 120min
Device Intensity	19 levels
Material	ABS
Stimulation modes	8
Size of Device (LxWxH)	Length & Width 49mm, Height 12.5mm

4. Instructions for operation and use

1. Connect the MASTOGO TENS Therapy device to the massage pads, peel off the protective film of the massage pads, and then paste the product on the part that needs to be massaged. Pay attention to keep the skin of the massage site clean and dry, and the massage pads should be evenly applied to the skin.
2. After short pressing the “M” button for 2 seconds to turn on the device, the default working time is 20 minutes, and the time cannot be adjusted. Press the “M” button to switch between modes. Please note that after switching modes, the intensity will automatically become 0 and the corresponding indicator light will be always on.
3. Press “+” or “-” to adjust the intensity. When working, the green light on the outer ring will be always on, the mode and the light will correspond one by one, and the intensity indicator light will show the current intensity. The green light is on, indicating that the device is enabled. (The current mode indicator flashes when there is no load, and it will return to normal display when there is a load.).
4. After the massage, short press the “M” button for 2 seconds to turn off the power (no load, no operation for 2 minutes, it will automatically shut down). Then slowly peel the edge of the massage pad from the skin. Finally, re-attach the protective film to the gel pad and store it.

5. Battery Charging, and Battery Status

1. Connect the USB cable to the TENS device.
2. Connect the USB cable with suitable adapters (not induced).
3. It takes 120 minutes to charge from low power to full power. The product cannot be turned on and there can be no output during charging. After charging, the indicator light flashes continuously between different modes.
4. When charging, the indicator light switches from mode 1 to mode 8, and all modes are on. When fully charged, unplug the charging cable and keep the device turned off.

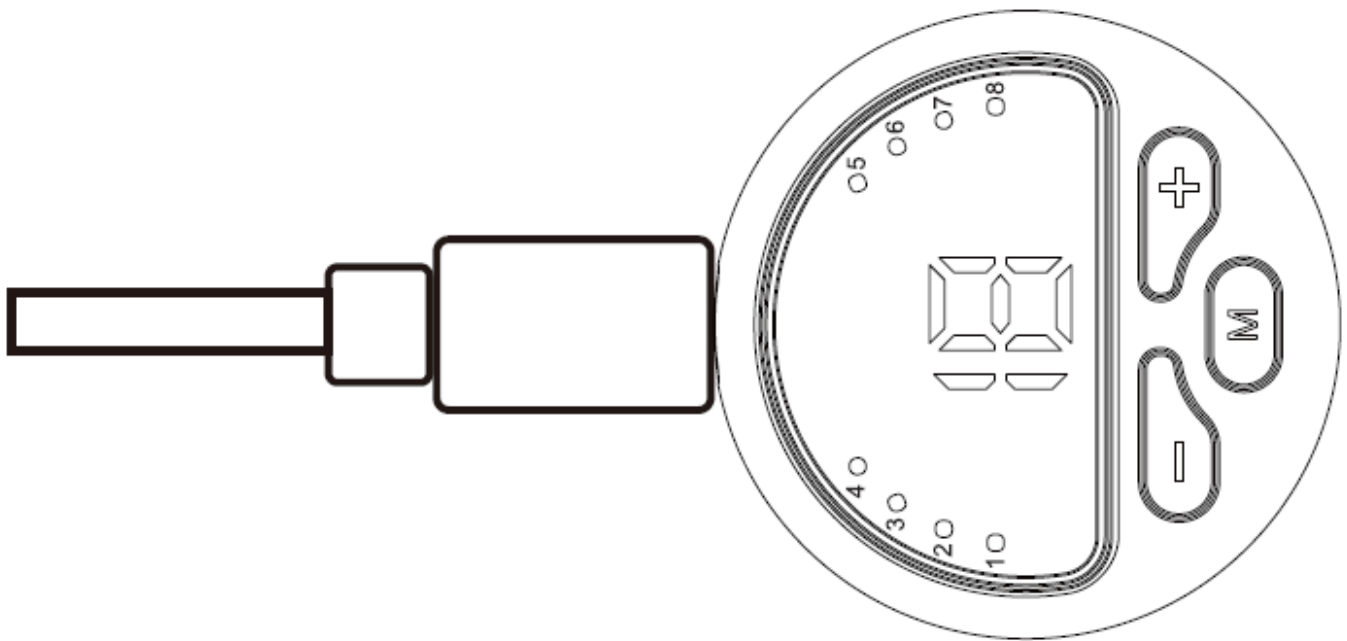6. Troubleshooting

Troubleshooting Guide

Problem	Cause	Solution
Can't feel any stimulation	The intensity level is too low.	Press the + button until you feel stimulation.
	The gel electrode patch is not firmly attached to the skin.	Press the gel electrode pad down against your skin to ensure even contact.
	The adhesive surface of the gel patch is damaged.	Replace it with a new gel pad.
	The adhesive surface of the gel patch is dirty or dusty.	See Maintenance and Disposal and clean the gel pad.
The electrical stimulation is too weak.	The battery is very low.	Charge the device.
	The adhesive surface of the gel patch is damaged.	Replace it with a new gel pad.
	The adhesive surface of the gel patch is dirty or dusty.	See Maintenance and Disposal and clean the gel pad.
The electrical stimulation is too strong.	The intensity level is too high or the gel patch is not properly attached.	Press the - button to decrease the intensity or ensure the gel pad is properly adhered.
The stimulation stopped unexpectedly.	The battery is low.	Charge the device.
The gel electrode patch is not adhering firmly to the skin.	The adhesive surface of the gel patch is dirty or dusty.	See Maintenance and Disposal and clean the gel pad.
	The adhesive surface of the gel patch is damaged.	Replace it with a new gel pad.

Irritated or red skin	The patch is not properly adhered.	Clean the skin and apply firmly (no creams or oils).
	The patch surface is dirty or dry.	Clean the patch (see Maintenance section).
	The gel patch is damaged.	Replace with a new gel patch.

7. Maintenance and disposal

7.1 Limited Warranty

This limited warranty guarantees the first purchaser of the TENS device the following: if the device fails to function within normal tolerances due to a defect in materials or workmanship within a period of 6 months from the date of delivery, you can contact us through the purchase channel or send us after-sales emails. We will resolve your issue promptly, so please rest assured.

7.2 Maintenance and Storage

If the device is dirty or has accumulated gel pad residue, use a damp cloth with mild soap to clean it.

Store the device and supplied items in their original packaging when not in use.

The device should be used, transported, and stored at temperatures between 10°C and 43°C (50°F and 110°F), with a relative humidity between 30% and 85%. All values are subject to a tolerance of +/- 10%.

The electrode gel pad should be replaced when it no longer adheres to the skin (it can be reused approximately 40-50 times). When you need a replacement, purchase electrode gel pads from our store.

Limitations of Use and Storage		
 500 hPa 1060 hPa Atmospheric pressure limitations	 50°F 110°F 10°C 43°C Temperature limitations	 30% 85% Humidity limitations

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