GADNIC

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MASAJ308

Índice

Español	5
Português	9
English	13

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SHIATSU MASSAGER WITH HEATER

**Instruction
manual**

Read these instructions carefully before use.


ATTENTION

1. **IMPORTANT:** Anyone who; may be pregnant or have recently given birth, has heart problems or a pacemaker, has diabetes or any other sensory problems, diabetes and/or thrombosis, arthritic dysfunction, injuries or inflammation, weak bones, osteoporosis, high fever, scoliosis, is being treated by a doctor, should consult a health professional before using this device.
2. Do not pull the lock or touch any of the internal parts for any reason.
3. Do not clean the switch with water.
4. When the product is in use or charging, make sure someone is present. When not in use or when leaving, unplug the power cable from the socket.
5. Do not use a damaged cable or plug.
6. Do not spray flammable oils, paints or pesticides.
7. Do not use in damp areas, such as bathrooms.
8. Avoid direct sunlight during long periods of use to prevent possible damage.
9. Do not open the structure or replace internal parts. Avoid short circuits and electric shocks.
10. Do not place near fire or combustibles.
11. Keep this product out of the reach of children.

Instructions for use

1. Insert the plug into the charging connector.
2. Pressing the power button  automatically starts the massage and heating function.
3. Press the button  to make the massage head move clockwise (or anti-clockwise).
4. This product has three levels of intensity. To adjust them, press the button  to set the speeds: low, medium and high.
5. Press the button  to desactivar the heater function; press again to activate it.
6. Press the power button  to switch off and remove the plug from the socket after using the equipment.

Schematic of the control panel and massage zones

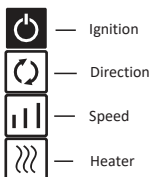Neck and shoulders

Back

Waist

Legs

NOTES

1. This product has a heating function, please pay attention to people who are sensitive to heat.
2. Children may only use this product under adult supervision.
3. If there is any discomfort or anomaly, discontinue use immediately.
4. After use, press the power switch to switch it off and unplug the power cable from the socket.
5. Drinking 500 ml of warm water or tea after the massage will help your body's metabolism.
6. If the cable is damaged, it must be replaced to avoid the risk of fire.

Storage

1. Unplug the equipment when not in use.
2. Do not spill water or any other liquid on the product.
3. Clean the product with a soft cloth.
4. If the product is dirty, clean it with a soft cloth dampened with water or a neutral detergent. Wipe off the dust with a dry cloth.
5. Do not clean the product with water.
6. Store the product in a cool, dry place. Do not expose the product to direct sunlight, high temperatures, fog or humid environments.
7. Do not use this product while using body sprays.

Common problems

Cause	Treatment
The device won't switch on	Confirm that the switch has been pressed.
The team is in the middle of a massage	1. Accidentally pressed the ignition switch. 2. The product has reached its usage limit of 15 minutes.

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