GADNIC

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KIT DE VENTOSAS

VENTO005


Índice

Español	5
Português	9
English	13

GADNIC

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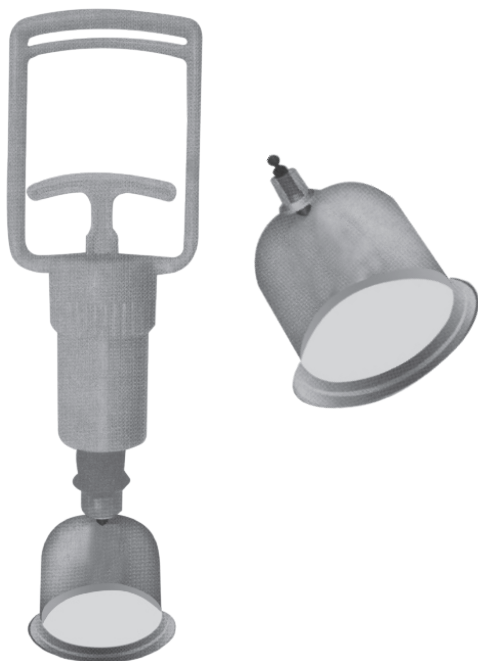1、 Descripción general

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
Flakes for vacuum

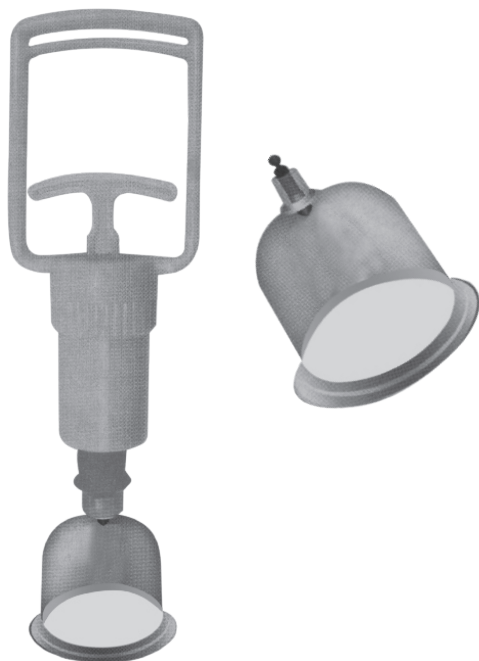1、General description

Our new suction cup instrument replaces the traditional technique of evacuating air by ignition with the vacuum air purge technique, which causes the suction cup instrument to be absorbed by the acupuncture points on the surface of the body to cause counterflow or slight sedimentation of the blood, in order to facilitate smooth flow in the channels and collaterals, promote blood circulation, eliminate diseases and thus make the qi and blood flourish. The suction cup instrument can produce a powerful magnetic field effect and is very effective in relieving pain, reducing swelling, lowering blood pressure, reducing fat in the blood, ending dysrhea, detoxifying, preventing the proliferation of cancer cells, as well as other efficacies such as eliminating accumulated fat and keeping fit. The suction cup instrument can be used for various types of medical treatment, such as internal medicine, surgery, gynaecology, paediatrics, treatment of skin diseases and diseases of the five sense organs. As a high-quality medical product, the new suction cup instrument has advantages such as a fine perspective, high safety and reliability, high efficiency, healing from the surface to the root, cost-effectiveness and durability.

2、Method of implementation


- (1) Pull the top bar of the cup valve to ensure a smooth flow of air before starting treatment;
- (2) Attach the vent nozzle to the top of the suction cup and pull the rod vertically about 4 times to create adequate negative pressure;
- (3) When the treatment is complete, gently pull up on the valve bar at the top of the suction cup to remove the instrument from the suction cup;
- (4) Use a disinfectant liquid to clean the suction cup instrument after application. Never immerse it in liquid.

3、Treatment methods

- (1) Choose suitable positions, such as lying on your stomach, on your stomach or on your side, depending on the state of your illness.
- (2) A treatment lasts 10-15 minutes, is carried out every day or every other day, and 7 treatments make up a treatment period. You can start the next treatment period after a 4-day break.
- (3) The suction cup tool can be used in various ways, for example using a single cup or several cups at the same time. In the case of the elderly or children, it is best to use fewer cups and pull the cup more gently, and the treatment time should be comparatively shorter, usually 5-10 minutes.
- (4) This product is especially effective in curing heatstroke (i.e. scraping therapy for heatstroke) and unfixed acupuncture points (i.e. any painful pressure point on the body).

4. Common cold, like hot


Symptom: Fever, distended headache, cough, thick phlegm, dry mouth, angina, pale skin and tongue. Selected acupuncture points: Dazhu, Quchi, Taiyang (fast suction cups), Fengchi, Tiantu. The effect will be best if the suction cup is moved over the bak du impulse and both urinary bladder channels until the skin turns red before placing the suction cup over the Yu acupoint and other relevant acupoints.

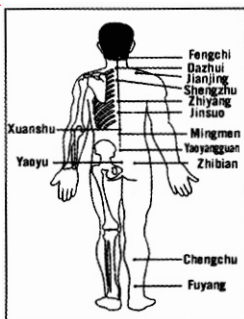5. Diseases of bone surgery

(1) Cervical spondylopathy

Symptom: Neck pain, coma, stiffness, limitation of movement or pain extending to the forearm, abductor or whirling or feeling unwell. Selected acupoints: Dazhui, Dazhu, Jianzhongyu, Jianwaiyu, Jianjing, Tiansong, Quchi, Waiguan and Fengchi.

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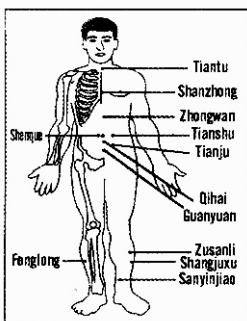
Symptom: Shoulder pain/pressure, limited movement, more intense pain at night. Acupuncture selected: Jianliao, Jianyu, Jianzhen, Tiansong, Dazhui, Fengchi, Quchi and Hegu.


(3) Chronic lumbago

Symptom: Continuous pain in the waist, the pain is aggravated by tiredness or cold. Selected acupuncture points: Shengyu, Yaoyangguan, Qihai, Mingmen, Huantiao, Weizhong, Yangling quan, Yongquan and other non-fixed acupuncture points.

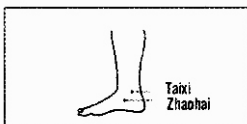
(4) Stiff collar

Symptom: neck stiffness or pain caused by incorrect sleeping posture or neck tension. Selected acupuncture points: Dazhui, Fengchi, Fengmen, Jianjing, Jianzhongyu, Jianwaiyu and other non-fixed acupuncture points.

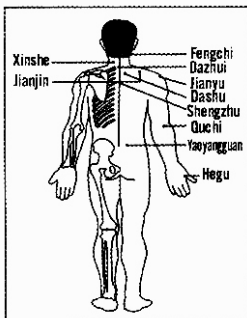(5) Soft tissue injuries

Symptom: Partial swelling, pain and limited movement caused by muscle ligament damage, etc. Upper limbs ; Selected acupuncture points: Hegu, Shousanli, Quchi and non-fixed acupuncture points. Selected acupuncture points: Yaoyangguan, Weizhong, Dazhui, Yongquan and other non-fixed acupuncture points. Selected acupuncture points for the lower limbs: Weizhong, Chengshan, Yanglingquan and other temporarily selected acupuncture points.

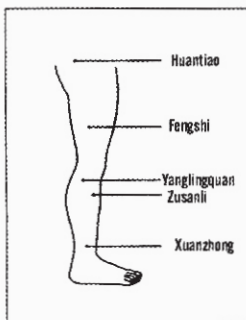
(6) Heel pain

Symptom: Pain in the heel, inability to stand or walk, caused by trauma, overexertion or heel spurs. Selected acupuncture points: Chengshan, Taixi, Kunlun, Yongquan, Zhaohai and other non-fixed acupuncture points.

(7) Bone nerve pain

Symptom: Throbbing, burning or widespread pain in the waist, hip, back of the thigh and calf or limitation of movement; the pain worsens after tiredness or cold. Selected points: waist reaction point, Shengyu, Yaoyangguan, Huantiao, Chibian, Jizhong, Chengshan, Yinmen, Xuanzhong, Fengshi, Yanlingquan.

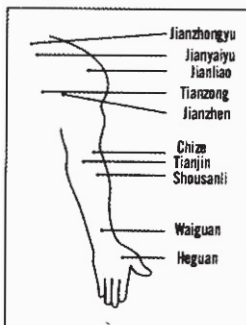6. Precautions

- (1) Treatment with suction cups should not be applied to severe heart disease, heart failure, diseases with a tendency to haemorrhage, patients with tumours, pregnant women, women with a lower abdomen in the menstrual period, elbow, trunk, heart position or in the event of a coma, bone fracture or seizure.
- (2) The room where the cupping treatment is carried out should be kept warm and avoid cold wind. It is advisable to drink lukewarm water after the cupping treatment and not to wash the cup marks on the body within 30 minutes of the cupping treatment.
- (3) The petechiae (blue, purple, black and red) that remain after the suction cups are placed are a type of symptomatic reaction. For example, black petechiae indicate cold, yellow and red petechiae indicate heat. Don't worry, though. The petechiae disappear 3 to 5 days later. It's normal to find blisters on the body for a certain period after treatment with suction cups. If the blister is large, disinfect it regularly and prick it with a disinfected needle, then clean the blister with a tampon and disinfect it with 75% medical alcohol to avoid infection.

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