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CEPILLO DENTAL
ELÉCTRICO

WATER010

INDICE

ESPAÑOL	03
PORTUGUÊS	06
ENGLISH	09

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
Nome de cada parte

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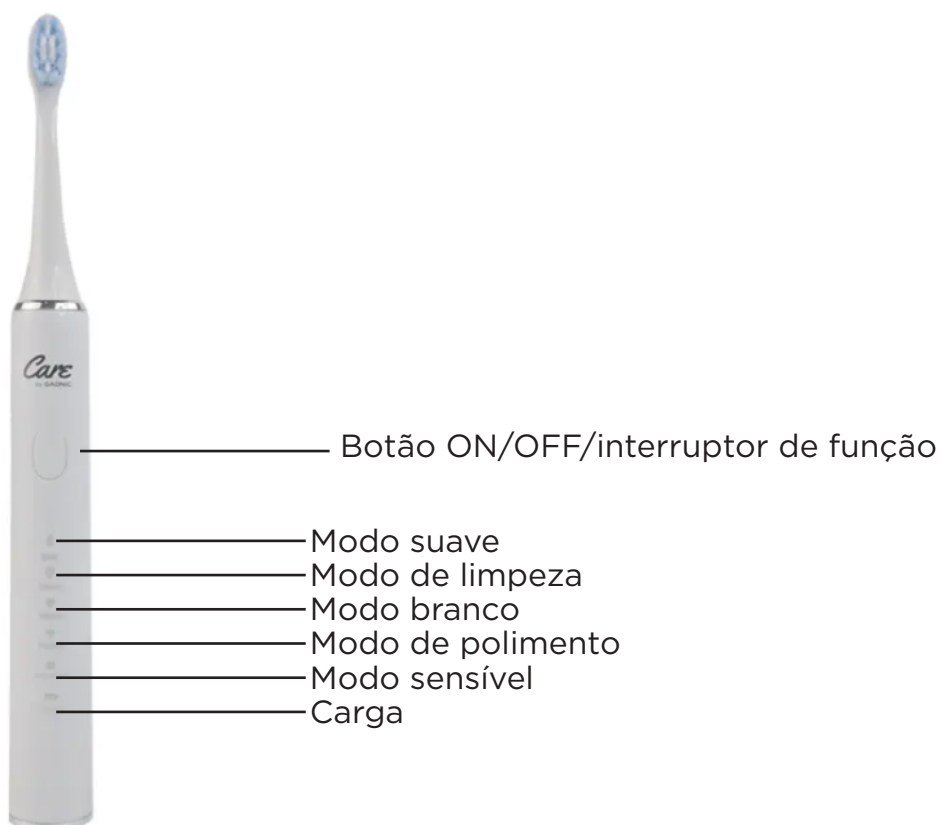This instruction manual contains safe operation and instructions for the correct use of this product. Read this manual carefully before use and take care not to damage or endanger the product.

Safety instructions

1. Please do not use it outside the specified use!
2. When connecting to the power supply, please check if the voltage is suitable for the local voltage.
3. Please use the original charger to avoid hazards.
4. Please avoid using abrasives with bleach and tooth stain remover together.
5. Please do not use with high intensity on teeth and gums.
6. Children and people who cannot operate by themselves or have increased oral sensitivity, please do not use.
7. If the teeth are being treated, or if the dental implant is not strong, do not use it.
8. Do not drop it, do not hit it violently.
9. Do not immerse in water for long periods of time.
10. This product is for personal use only and is not available for dental clinics or multi-patient organizations.

Warning

Children and persons with physical disabilities should use under the supervision of a guardian.

In case of periodontal disease and oral symptoms, use under the supervision of a physician.

If you feel discomfort or pain in your mouth, stop using it and consult your dentist. In case of electric toothbrush failure, consult the local administration on how to dispose of the electric toothbrush to protect the environment as much as possible.

CAUTION! BEFORE THE FIRST USE YOU MUST CHARGE THE PRODUCT FOR 6 HOURS. DURING THE CHARGING PROCESS A GREEN LIGHT WILL BE LIT, HOWEVER, COMPLETE THE INDICATED CHARGING TIME.

How to use?

1. Moisten the toothbrush head and apply an adequate amount of toothpaste.
2. Place the bristles on the teeth at a 45-degree angle aligned with the gum line.
3. Press the on/off button, the brush will start or stop working.
4. You just need to gently move the bristles to the place you need to clean. The pressure should be less than 150g (Sanie toothbrush relies on the high frequency vibration of the bristles to clean the mouth, so when brushing

the teeth, the using force is much less than that of a normal toothbrush, about 1/8).

5. Dentists advise dividing the oral cavity into 4 areas (as shown above), brushing slowly to ensure that each area has approximately 30 seconds of cleaning time. The machine has a 30 second area reminder function, it will pause for 0.5 seconds every 30 seconds to remind you to change areas to brush your teeth.

6. Press the on/off button, the brush will start working, then press the button for three seconds to stop working. This toothbrush has five brushing modes: Gentle, Clean, White, Polished, Sensitive. You can press to select the working mode according to your needs. This toothbrush also has a working mode memory function, the next time you use the brush, it will directly enter the mode you used, no need to adjust.

Product features

1. High frequency sonic vibration: 33,000 times per minute high frequency sonic vibration, in only 15 days can remove stains (including tea stains, smoke stains), achieving whitening effect.


2. Ergonomic design: waterproof, non-slip, comfortable feeling, safe and reliable.

3. USB fast charging design: Unplug the waterproof plug on the bottom of the toothbrush. Plug the DC connector of the charging cable into the charging port on the bottom of the toothbrush before inserting, and then insert the USB connector into the computer or 5V USB socket. The maximum charging time is 6 hours. The green light flashes slowly when it is charging, and the green light turns on when it is fully charged.

4. Three brush heads: comes with three soft brush heads, which can be replaced anytime, anywhere, for long-term care of your oral health.

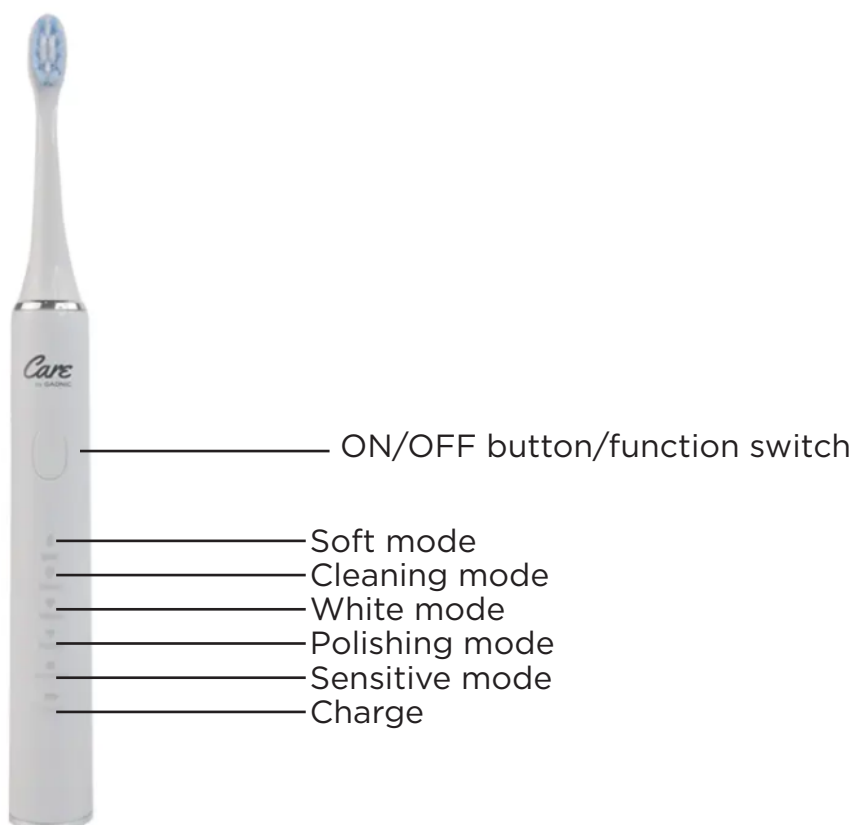5. Healthier brush safety cap: set the brush safety cap, not only to prevent bacteria invasion, but also more convenient and healthier to carry.

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