GADNIC

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MINI MICROSCOPIO PORTATIL

MICROS25

MICROSCOPIO

PORTATIL

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
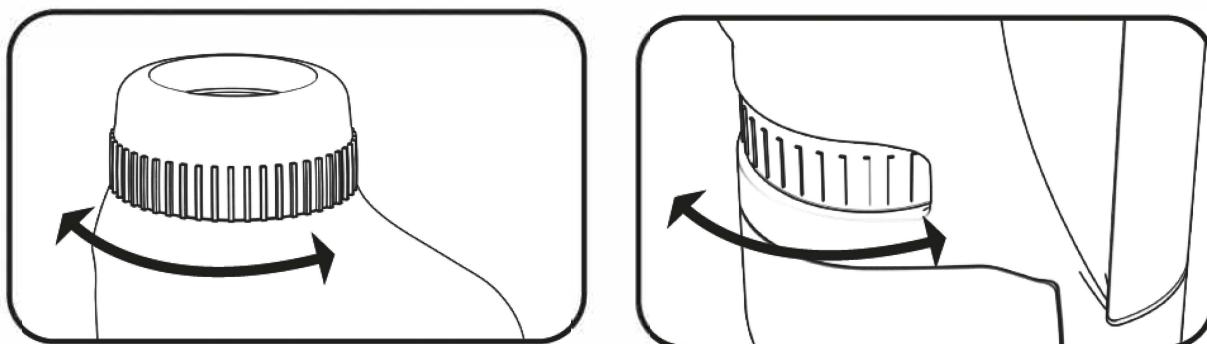
PORTABLE

User's manual

PLEASE READ THIS MANUAL CAREFULLY
BEFORE USING THE PRODUCT.

Within a small world

The world in which all living things live is very mysterious and vast. People's eyes can see many beautiful things in the world, for example, beautiful trees, flowers and plants, imposing buildings, trains and cars; but there are many hidden objects in the world that people cannot see because they are too small. Our eyes cannot see the so-called microscopic world.

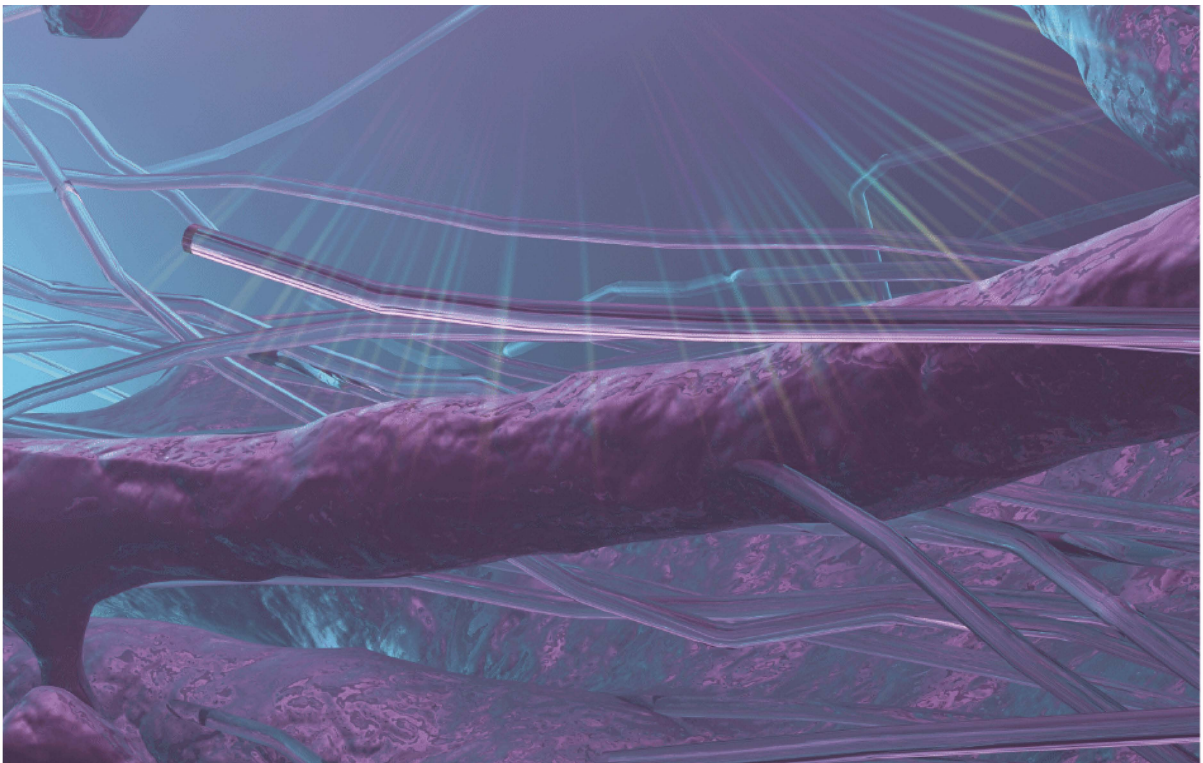Better observation tools

How can we observe the world in more detail? This is not difficult to do at home. You can use the magnifying glass to observe something that cannot be seen clearly with the naked eye, such as using it to observe insects, for example.

So how to observe objects smaller than insects? Scientists combined two magnifying glasses to make a microscope, and discovered that many things that were previously invisible to the naked eye could be observed under it: small insects swimming in water, stomata on leaves, and so on.

So how to observe objects smaller than insects? Scientists combined two magnifying glasses to make a microscope, and discovered that many things that were previously invisible to the naked eye could be observed under it: small insects swimming in water, stomata on leaves, and so on.

The world under the microscope

What does the world look like under the microscope?
Let's check it out!


Onion epidermal cells


Apple epidermal cells


Plant cells

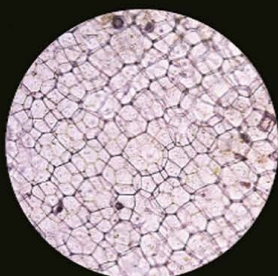
Brahman ginseng
cell

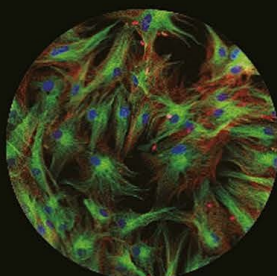
Skin cell

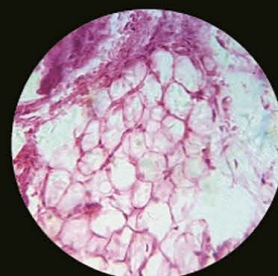Skeletal cells

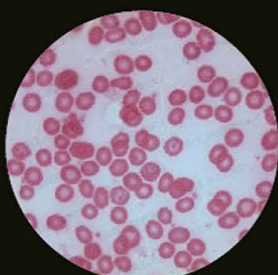
Human red
blood cells

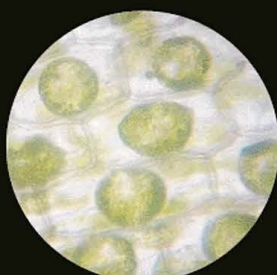Plant cell


Purple cell

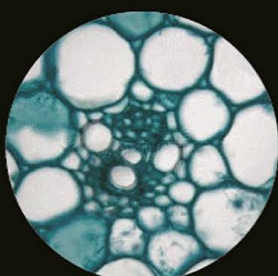Corn cell

Histocyte

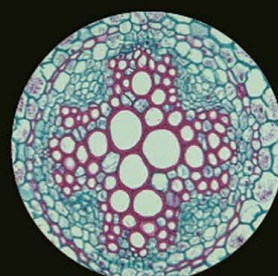Leaf cell

Operating Instructions

Simple version

(Product requires 2 "AAA" batteries)

PRECAUTIONS FOR THE USE OF BATTERIES:


1. ALKALINE, CARBON-ZINC AND OTHER TYPES OF BATTERIES CAN BE USED. DO NOT MIX ALKALINE, STANDARD (CARBON-ZINC) OR RECHARGEABLE (NICKEL-CADMIUM) BATTERIES.
2. NON-RECHARGEABLE BATTERIES CANNOT BE CHARGED.
3. RECHARGEABLE BATTERIES MAY ONLY BE CHARGED UNDER ADULT SUPERVISION.
4. RECHARGEABLE BATTERIES MUST BE REMOVED FROM THE TOY BEFORE CHARGING.
5. DO NOT MIX DIFFERENT TYPES OF BATTERIES OR OLD AND NEW BATTERIES.
6. BATTERIES MUST BE INSERTED WITH THE CORRECT POLARITY (+, -).
7. EXHAUSTED BATTERIES MUST BE REMOVED FROM THE TOY.
8. PLEASE REMOVE THE BATTERIES WHEN THE PRODUCT IS NOT USED FOR A LONG TIME.
9. DO NOT THROW THE BATTERIES INTO FIRE.

WARNING

ALL PACKAGING MATERIALS, SUCH AS TAPE, PLASTIC SHEETING, STRING, LABELS, ETC., ARE NOT PART OF THE PRODUCT AND ARE NOT SUITABLE FOR CHILDREN TO PLAY WITH. FOR THE SAFETY OF YOUR CHILDREN, HANDLE THE ABOVE MENTIONED PACKAGING MATERIALS PROPERLY. THE USER'S MANUAL CONTAINS IMPORTANT INFORMATION. PLEASE KEEP IT FOR FUTURE USE.

INSTALLATION OF THE BATTERIES:


1. CONFIRM THAT THE PRODUCT IS TURNED OFF.
2. PRESS THE LEFT AND RIGHT SIDES OF THE BATTERY COVER HOLE OF THE PRODUCT AND THEN OPEN THE BATTERY COVER BACKWARDS.
3. INSTALL 2 NEW 1.5V AAA BATTERIES WITH THE CORRECT POLARITY (+, -) AS SHOWN IN THE FIGURE ON THE RIGHT (NON-RECHARGEABLE BATTERIES).
4. CLOSE THE BATTERY DRAWER COVER. THEN OPEN THE BATTERY COVER BACKWARDS.
5. INSTALL 2 NEW 1.5V AAA BATTERIES WITH THE CORRECT POLARITY (+, -), AS SHOWN IN THE FIGURE ON THE RIGHT (NON-RECHARGEABLE BATTERIES).
6. CLOSE THE BATTERY DRAWER COVER.


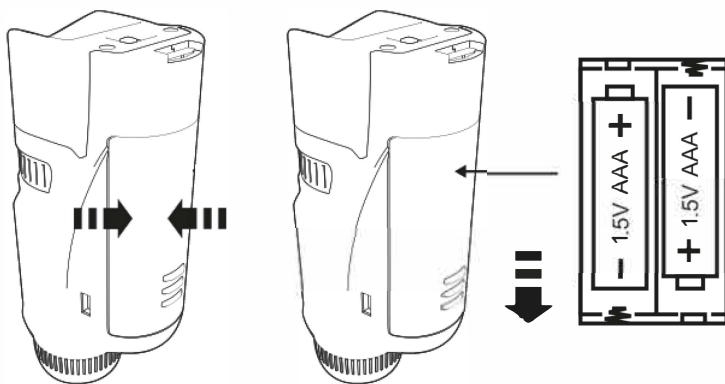
INSTRUCTIONS FOR USE

Manual observation mode

Step 1: Turn on the led light, point your eyes towards the ocular observation port and point the objective lens towards the observation object.

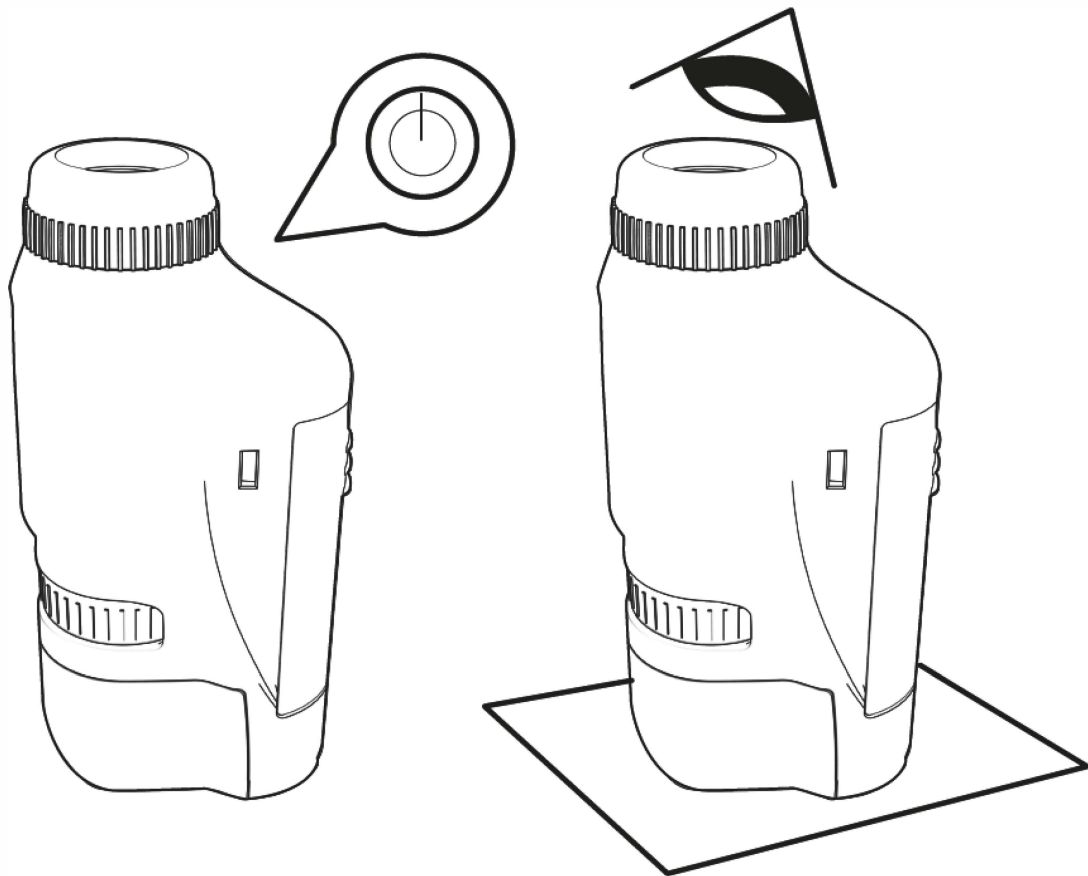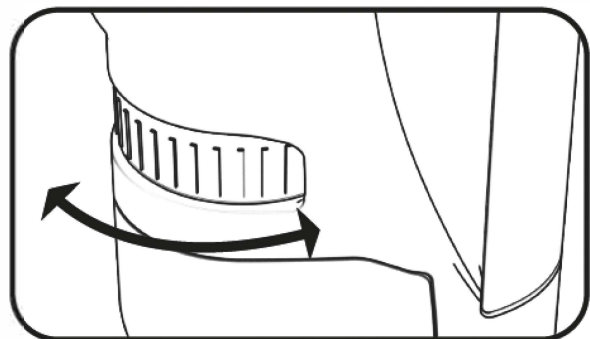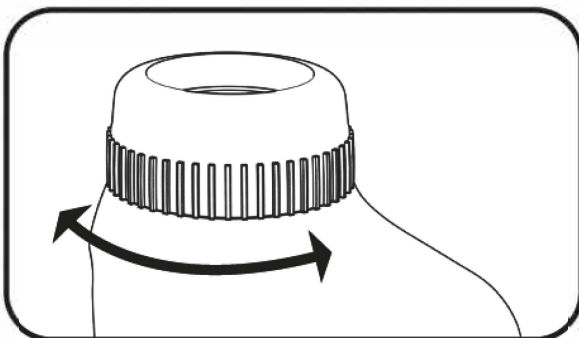Step 2: Manually adjust the zoom wheel until the image is clear.

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